

HPLC Application

ID No.: 18850

Pain panel extracted from human urine at cutoff - Kinetex C18 2.6µm 50x2.1mm LCMSMS

Column: Kinetex[®] 2.6 µm C18 100 Å, LC Column 50 x 2.1 mm, Ea

Dimensions: 50 x 2.1 mm ID

Order No: 00B-4462-AN

Elution Type: Gradient

Eluent A: 0.1% Formic acid in H₂O

Eluent B: 0.1% Formic acid in ACN

Gradient Profile:	Step No.	Time (min)	Pct A	Pct B
	1	0	95	5
	2	1.3	40	60
	3	1.31	5	95
	4	1.51	5	95
	5	1.52	95	5
	6	2.2	95	5


Flow Rate: 1.2 mL/min

Col. Temp.: ambient

Detection: Mass Spectrometer (MS) @ amu (ambient)



Products used in this application:

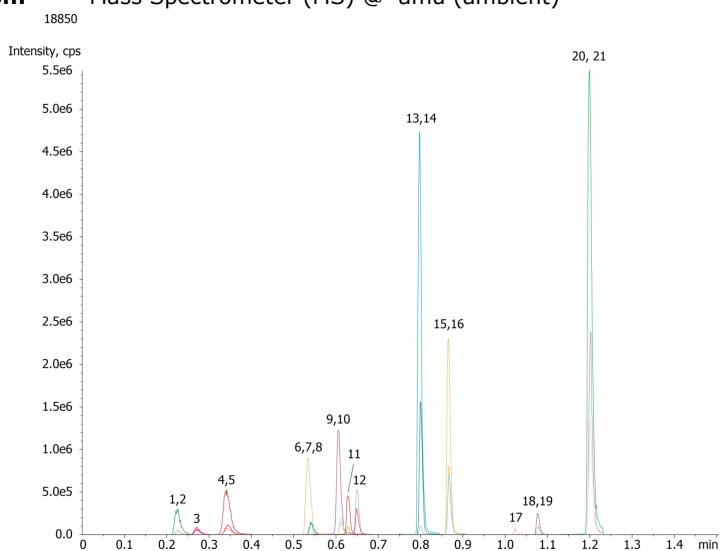

Pain panel extracted from human urine at cutoff - Kinetex C18 2.6µm 50x2.1mm LCMSMS

ANALYTES:

- 1 Morphine
- 2 Morphine-d3
- 3 Oxymorphone
- 4 Hydromorphone
- 5 Hydromorphone-d6
- 6 Codeine
- 7 Codeine-d6
- 8 Naloxone
- 9 Oxycodone
- 10 Oxycodone-d6
- 11 Naltrexone
- 12 Hydrocodone
- 13 cis-Tramadol
- 14 C13-D3-tramadol
- 15 Meperidine
- 16 Meperidine-d4
- 17 Fentanyl
- 18 Buprenorphine
- 19 Buprenorphine-d4
- 20 Methadone
- 21 Methadone-d3

