

# HPLC Application

ID No.: 18866

## Terfenadine on Lux 5 µm Cellulose-1 by Chiral SFC (loading and scale up)

**Column:** Lux<sup>®</sup> 5 µm Cellulose-1, LC Column 250 x 4.6 mm, Ea

**Dimensions:** 250 x 4.6 mm ID

**Order No:** 00G-4459-E0

**Elution Type:** Isocratic

**Eluent A:** Methanol with 0.1% DEA

**Eluent B:** CO<sub>2</sub>

Gradient Profile:	Step No.	Time (min)	Pct A	Pct B
	1	0	25	75

**Flow Rate:** 2.5 mL/min

**Col. Temp.:** 35 °C

**Detection:** UV-Vis Abs.-Variable Wave.(UV) @ 220 nm (ambient)


**Analyst Note:** Application using combination of UV-Vis detector (upper trace) and PDR-Chiral Advanced Laser Polarimeter (lower trace) © 2024 Phenomenex, Inc.

1.5mg in 50µL on-column injection of Terfenadine.

18866



Products used in this application:

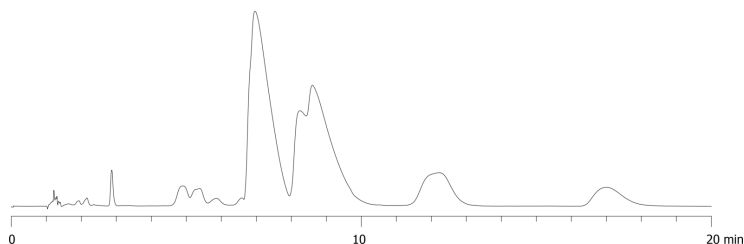

### ANALYTES:

1 Terfenadine

2 Terfenadine

