

# HPLC Application

ID No.: 18932

## Hydrochlorothiazide Related Substance Analysis on Luna 5u C18 250 x 4.6mm ID

**Column:** Luna® 5 µm C18(2) 100 Å, LC Column 250 x 4.6 mm, Ea

**Dimensions:** 250 x 4.6 mm ID

**Order No:** 00G-4252-E0

**Elution Type:** Isocratic

**Eluent A:** 100mM Na Phosphate Buffer pH=3.0

**Eluent B:** Acetonitrile

Gradient Profile:	Step No.	Time (min)	Pct A	Pct B
	1	0	90	10

**Flow Rate:** 1.5 mL/min

**Col. Temp.:** ambient

**Detection:** UV-Vis Abs.-Variable Wave.(UV) @ 254 nm (ambient)

**Detector Info:** MBP=190 Bar


**Analyst Note:** SecurityGuard™ Guard Cartridge System extends column lifetime.

- SecurityGuard Cartridges, C18 4 x 3.0mm ID, 10/Pk Part No.: AJ0-4287

- Holder Part No.: KJ0-4282



Products used in this application:

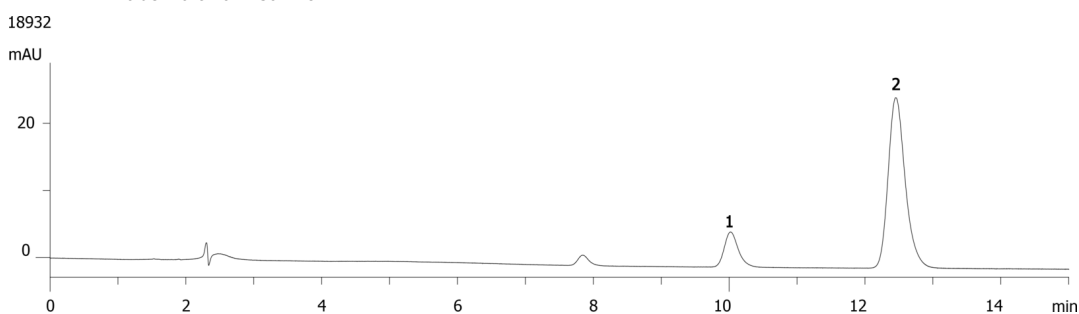

### ANALYTES:

- 1** Chlorothiazide  
Retention Time: 10.02 min
- 2** Hydrochlorothiazide  
Retention Time: 12.45 min

