


Impact of column length on Chiral Separations: 250mm length

Column: N/A
Dimensions: N/A
Order No: N/A
Elution Type: N/A
Eluent A: Ethanol with 0.1% DEA
Eluent B: Carbon Dioxide
Gradient

| Step No. | Time (min) | Pct A | Pct B |
|----------|------------|-------|-------|
| 1 | 0 | 25 | 75 |

Profile:
Flow Rate: 2.5 mL/min
Col. Temp.: 35 °C
Detection: UV-Vis Abs.-Fixed Wave. (UV) @ 254 nm (35 °C)
Detector Info: PDA

Analyst Note: SecurityGuard™ Guard Cartridge System extends column lifetime.
- SecurityGuard Cartridges, Silica 4 x 3.0mm, 10/Pk Part No.: AJ0-4348
- Holder Part No.: KJ0-4282

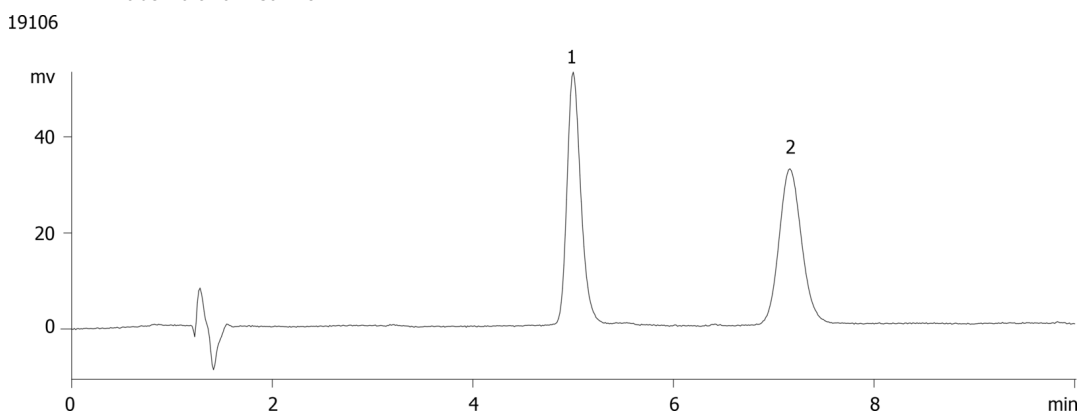

ANALYTES:

1 Propafenone

