


## Impact of column length on Chiral Separations: 250mm length

**Column:** N/A  
**Dimensions:** N/A  
**Order No:** N/A  
**Elution Type:** N/A  
**Eluent A:** Ethanol with 0.1% DEA  
**Eluent B:** Carbon Dioxide  
**Gradient**

Step No.	Time (min)	Pct A	Pct B
1	0	25	75

**Profile:**  
**Flow Rate:** 2.5 mL/min  
**Col. Temp.:** 35 °C  
**Detection:** UV-Vis Abs.-Fixed Wave. (UV) @ 254 nm (35 °C)  
**Detector Info:** PDA  
**Analyst Note:** SecurityGuard™ Guard Cartridge System extends column lifetime.

- SecurityGuard Cartridges, Silica 4 x 3.0mm, 10/Pk Part No.: AJ0-4348
- Holder Part No.: KJ0-4282

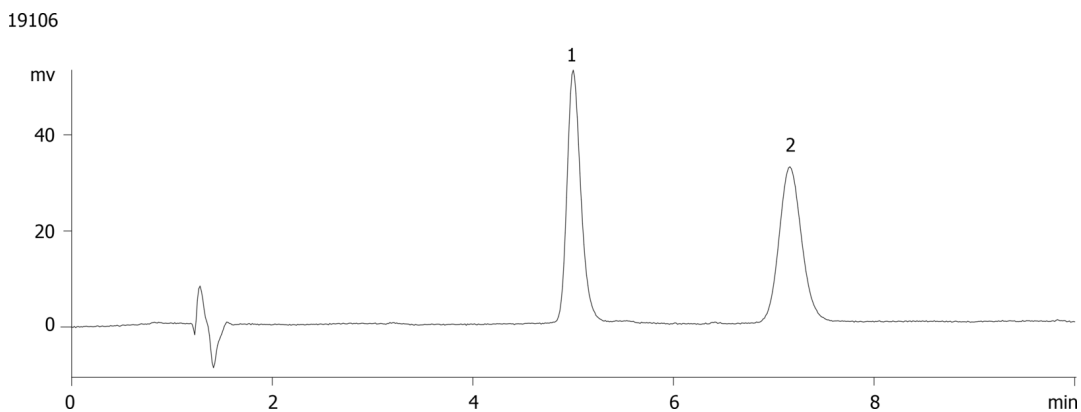

### ANALYTES:

- 1 Propafenone

