

USP Norethindrone, Progesterone and Mestranol on Luna 3 μ m C8(2) 150x4.60mm

Column: Luna® 3 μ m C8(2) 100 Å, LC Column 150 x 4.6 mm, Ea

Dimensions: 150 x 4.6 mm ID

Order No: 00F-4248-E0

Elution Type: Isocratic

Eluent A: 50:50 Acetonitrile:Water

Gradient	Step No.	Time (min)	Pct A
Profile:	1	30	100

Flow Rate: 1 mL/min

Col. Temp.: ambient

Detection: UV-Vis Abs.-Variable Wave.(UV) @ 200 nm (ambient)


Analyst Note: SecurityGuard™ Guard Cartridge System extends column lifetime.

- SecurityGuard Cartridges, C8 4 x 3.0mm ID, 10/Pk Part No.: AJ0-4290

- Holder Part No.: KJ0-4282



Products used in this application:

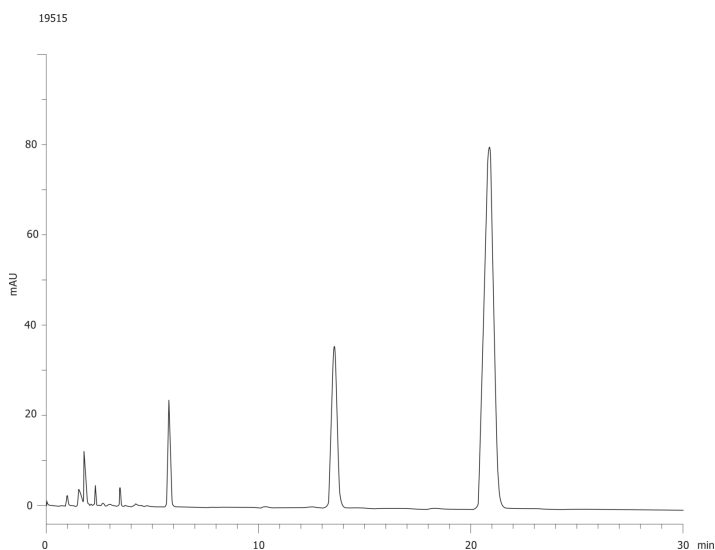

ANALYTES:

- 1** Norethindrone
Retention Time: 5.799 min
- 2** Progesterone
Retention Time: 13.584 min
- 3** Mestranol
Retention Time: 20.871 min

