

HPLC Application

ID No.: 19814

Degraded Myoglobin on Aeris WIDEPORE 3.6u XB-C18

Column: Aeris[™] 3.6 µm WIDEPORE XB-C18 200 Å, LC Column 150 x 2.1 mm, Ea

Dimensions: 150 x 2.1 mm ID

Order No: 00F-4482-AN

Elution Type: Gradient

Eluent A: 0.1% TFA in Water

Eluent B: Acetonitrile

Gradient Profile:	Step No.	Time (min)	Pct A	Pct B
	1	0	90	10
	2	5	90	10
	3	55	50	50

Flow Rate: 0.2 mL/min

Col. Temp.: ambient


Detection: UV-Vis Abs.-Variable Wave.(UV) @ 214 nm (ambient)

Analyst Note: SecurityGuard[™] Guard Cartridge System extends column lifetime.

- SecurityGuard ULTRA Cartridges, UHPLC WIDEPORE C18 for 2.1mm ID, 3/Pk Part No.: AJ0-8783

- Holder Part No.: AJ0-9000

AppID 19814

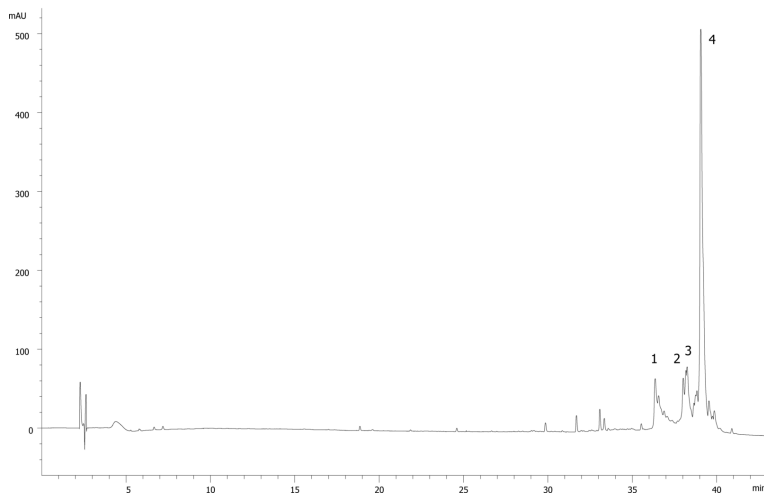

ANALYTES:

1 Myoglobin

Retention Time: 39 min



Products used in this application:

