

# HPLC Application

ID No.: 19822

## Degraded Myoglobin on Aeris WIDEPORE 3.6u XB-C18 150x2.1 Zoom In

**Column:** Aeris<sup>™</sup> 3.6  $\mu$ m WIDEPORE XB-C18 200  $\text{\AA}$ , LC Column 150 x 2.1 mm, Ea

**Dimensions:** 150 x 2.1 mm ID

**Order No:** 00F-4482-AN

**Elution Type:** Gradient

**Eluent A:** 0.1% TFA

**Eluent B:** Acetonitrile

Gradient Profile:	Step No.	Time (min)	Pct A	Pct B
	1	0	90	10
	2	5	90	10
	3	50	50	50

**Flow Rate:** 0.2 mL/min

**Col. Temp.:** ambient


**Detection:** UV-Vis Abs.-Diode Array (PDA) @ (ambient)

**Analyst Note:** SecurityGuard<sup>™</sup> Guard Cartridge System extends column lifetime.

- SecurityGuard ULTRA Cartridges, UHPLC WIDEPORE C18 for 2.1mm ID, 3/Pk Part No.: AJ0-8783

- Holder Part No.: AJ0-9000

AppID 19822

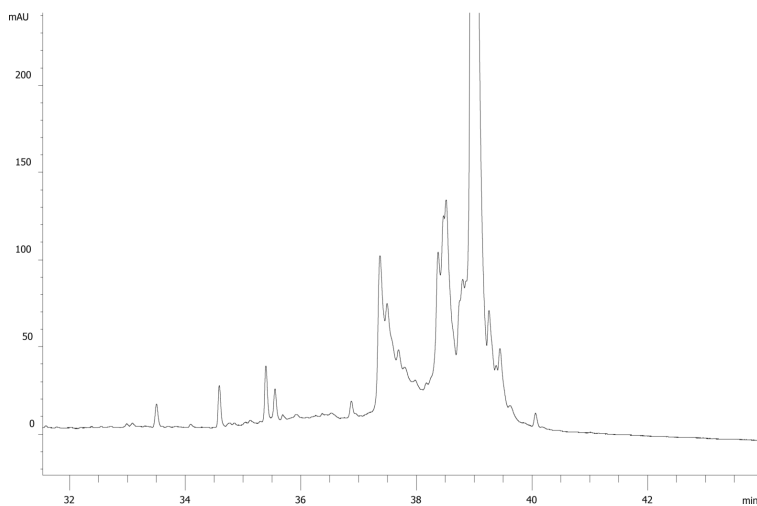

### ANALYTES:

1 Myoglobin



Products used in this application:

