

HPLC Application

ID No.: 19824

Degraded Myoglobin on Aeris WIDEPORE 3.6u XB-C18 250x2.1 Zoom In

Column: Aeris™ 3.6 µm WIDEPORE XB-C18 200 Å, LC Column 250 x 2.1 mm, Ea

Dimensions: 250 x 2.1 mm ID

Order No: 00G-4482-AN

Elution Type: Gradient

Eluent A: 0.1% TFA

Eluent B: Acetonitrile

Gradient Profile:	Step No.	Time (min)	Pct A	Pct B
	1	0	90	10
	2	5	90	10
	3	50	50	50

Flow Rate: 0.2 mL/min


Col. Temp.: ambient

Detection: UV-Vis Abs.-Diode Array (PDA) @ (ambient)


Analyst Note: SecurityGuard™ Guard Cartridge System extends column lifetime.

- SecurityGuard ULTRA Cartridges, UHPLC WIDEPORE C18 for 2.1mm ID, 3/Pk Part No.: AJ0-8783

- Holder Part No.: AJ0-9000

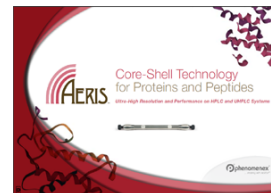

Products used in this application:

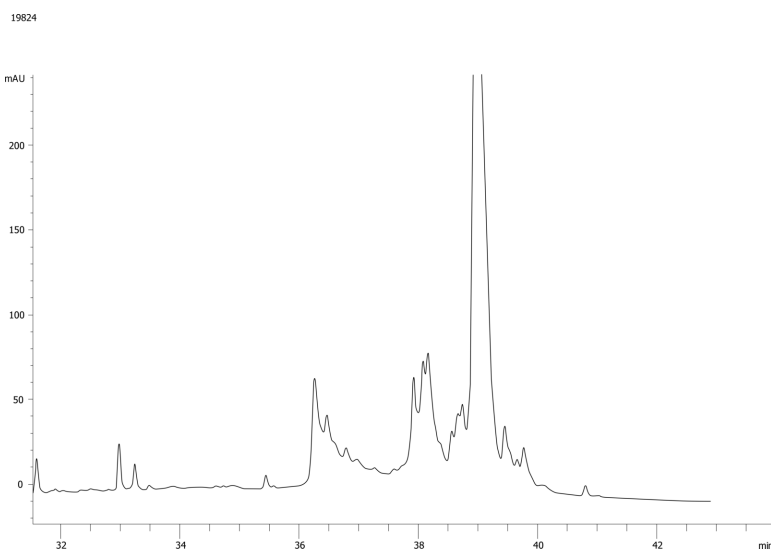

ANALYTES:

1 Myoglobin

