

# HPLC Application

ID No.: 19825

## Degraded Myoglobin on Aeris WIDEPORE 3.6u XB-C18 250x2.1 Zoom In High Flow

**Column:** Aeris™ 3.6 µm WIDEPORE XB-C18 200 Å, LC Column 250 x 2.1 mm, Ea

**Dimensions:** 250 x 2.1 mm ID

**Order No:** 00G-4482-AN

**Elution Type:** Gradient

**Eluent A:** 0.1% TFA

**Eluent B:** Acetonitrile

Gradient Profile:	Step No.	Time (min)	Pct A	Pct B
	1	0	90	10
	2	5	90	10
	3	50	50	50

**Flow Rate:** 0.3 mL/min

**Col. Temp.:** ambient


**Detection:** UV-Vis Abs.-Diode Array (PDA) @ (ambient)

**Analyst Note:** SecurityGuard™ Guard Cartridge System extends column lifetime.

- SecurityGuard ULTRA Cartridges, UHPLC WIDEPORE C18 for 2.1mm ID, 3/Pk Part No.: AJ0-8783


- Holder Part No.: AJ0-9000

AppID: 19825

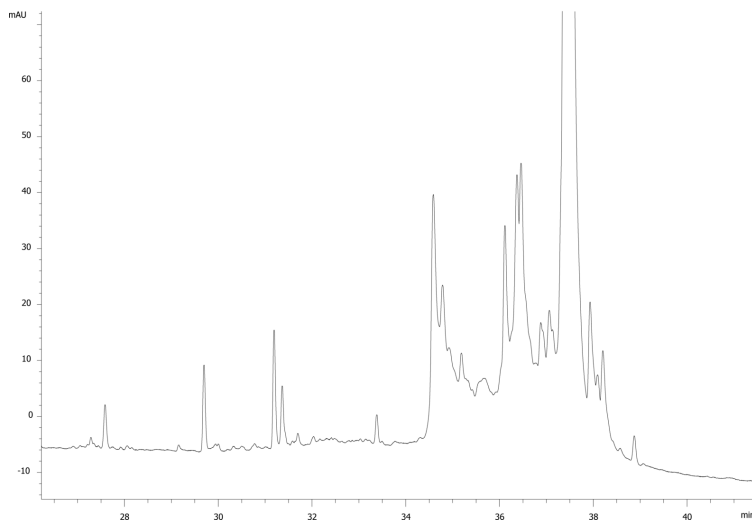

### ANALYTES:

1 Myoglobin

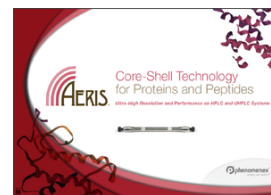

Products used in this application:

