

HPLC Application

ID No.: 19834

Myoglobin Degradation Aeris 3.6u WIDEPORE 150x2.1mm

Column: Aeris[™] 3.6 μ m WIDEPORE XB-C18 200 \AA , LC Column 150 x 2.1 mm, Ea

Dimensions: 150 x 2.1 mm ID

Order No: 00F-4482-AN

Elution Type: Gradient

Eluent A: Water with 0.1% TFA

Eluent B: Acetonitrile with 0.1% TFA

Gradient Profile:	Step No.	Time (min)	Pct A	Pct B
	1	0	90	10
	2	3	90	10
	3	33	35	65

Flow Rate: 0.3 mL/min

Col. Temp.: ambient

Detection: UV-Vis Abs.-Variable Wave.(UV) @ nm (ambient)


Analyst Note: SecurityGuard[™] Guard Cartridge System extends column lifetime.

- SecurityGuard ULTRA Cartridges, UHPLC WIDEPORE C18 for 2.1mm ID, 3/Pk Part No.: AJ0-8783

- Holder Part No.: AJ0-9000



Products used in this application:

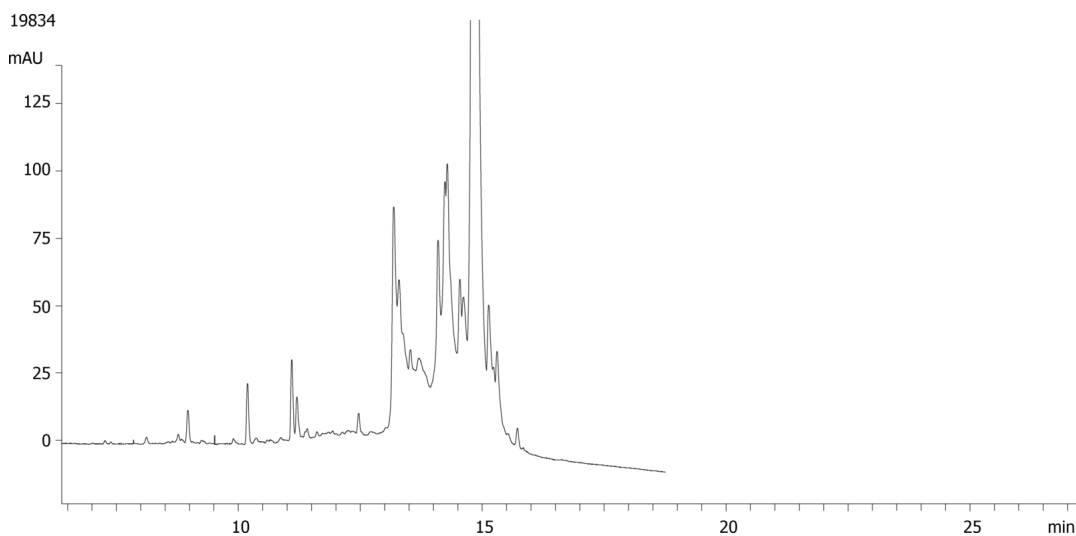

ANALYTES:

- 1 Myoglobin degradants

