

HPLC Application

ID No.: 19836

Myoglobin degradation: Aeris WIDEPORE 3.6u XB-C18 150x4.6mm (vs19837)

Column: Aeris[™] 3.6 μ m WIDEPORE XB-C18 200 \AA , LC Column 150 x 4.6 mm, Ea

Dimensions: 150 x 4.6 mm ID

Order No: 00F-4482-E0

Elution Type: Gradient

Eluent A: Water with 0.1% TFA

Eluent B: Acetonitrile with 0.1% TFA

Gradient Profile:	Step No.	Time (min)	Pct A	Pct B
	1	0	97	3
	2	3	97	3
	3	33	35	65

Flow Rate: 1.5 mL/min


Col. Temp.: 40 °C

Detection: UV-Vis Abs.-Variable Wave.(UV) @ nm (ambient)

Analyst Note: SecurityGuard[™] Guard Cartridge System extends column lifetime.

- SecurityGuard ULTRA Cartridges, UHPLC WIDEPORE C18 for 4.6mm ID, 3/Pk Part No.: AJ0-8769


- Holder Part No.: AJ0-9000

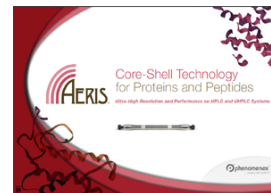

Products used in this application:



19836

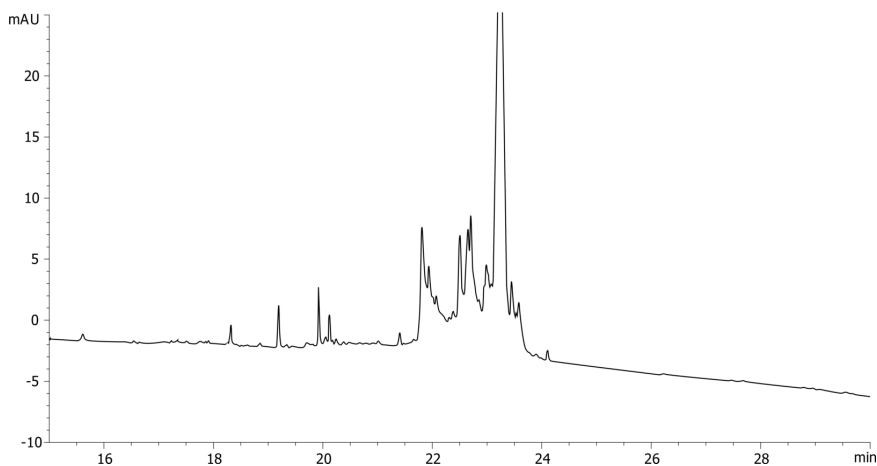

ANALYTES:

- 1 Myoglobin degradants

