

HPLC Application

ID No.: 19941

Pharmaceutical Drug Mix on Kinetex C18 2.6u 50 x 2.1mm @ 1.4 mL/min(Quotien Data)

Column: Kinetex® 2.6 µm C18 100 Å, LC Column 50 x 2.1 mm, Ea

Dimensions: 50 x 2.1 mm ID

Order No: 00B-4462-AN

Elution Type: Gradient

Eluent A: 0.1% formic acid in acetonitrile

Eluent B: 0.1% formic acid

Gradient Profile:	Step No.	Time (min)	Pct A	Pct B
	3	1.8	90	10
	4	1.9	90	10
	5	1.91	5	95
	6	2.5	5	95

Flow Rate: 1.4 mL/min




Col. Temp.: ambient

Detection: Mass Spectrometer (MS) @ amu (ambient)

19941



Products used in this application:

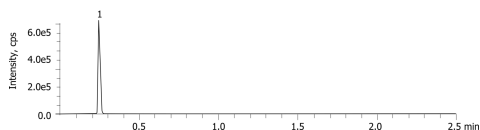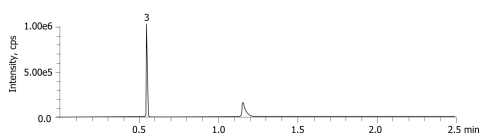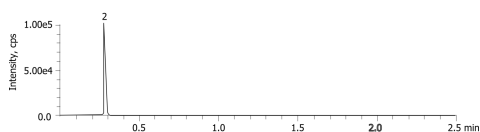

ANALYTES:

1 Cimetidine

Retention Time: 0.243 min

2 Acetaminophen

Retention Time: 0.278 min

3 Pindolol

Retention Time: 0.545 min

