

# HPLC Application

ID No.: 19946

## Pharmaceutical Drug Mix on Kinetex 2.6 µm C18 50 x 2.1mm @ 0.8 mL/min(Quotient Data)

**Column:** Kinetex® 2.6 µm C18 100 Å, LC Column 50 x 2.1 mm, Ea

**Dimensions:** 50 x 2.1 mm ID

**Order No:** 00B-4462-AN

**Elution Type:** Gradient

**Eluent A:** 0.1% formic acid in acetonitrile

**Eluent B:** 0.1% formic acid

Gradient Profile:	Step No.	Time (min)	Pct A	Pct B
	1	0	5	95
	2	0.25	5	95
	3	1.8	90	10
	4	1.9	90	10
	5	1.91	5	95
	6	2.5	5	95

**Flow Rate:** 0.8 mL/min




**Col. Temp.:** ambient

**Detection:** Mass Spectrometer (MS) @ amu (ambient)

19946



Products used in this application:

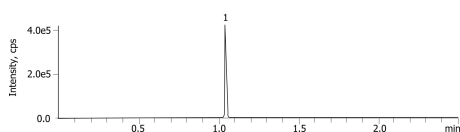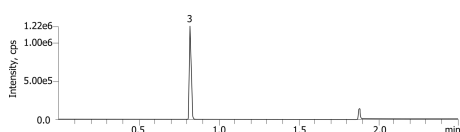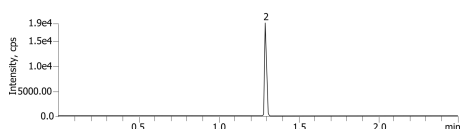

### ANALYTES:

**1** Diltiazem

Retention Time: 1.04 min

**2** Terfenadine

Retention Time: 1.29 min

**3** Acebutolol

Retention Time: 0.824 min

