

# HPLC Application

ID No.: 19978

## Isocratic Chromatogram of Kinetex without SecurityGuard ULTRA

**Column:** Kinetex® 2.6 µm XB-C18 100 Å, LC Column 50 x 2.1 mm, Ea

**Dimensions:** 50 x 2.1 mm ID

**Order No:** 00B-4496-AN

**Elution Type:** Isocratic

**Eluent A:** 50

**Eluent B:** 50

Gradient Profile:	Step No.	Time (min)	Pct A	Pct B
	1	0	100	0


**Flow Rate:** 0.5 mL/min

**Col. Temp.:** ambient

**Detection:** UV-Vis Abs.-Variable Wave.(UV) @ nm (ambient)  
19978



Products used in this application:

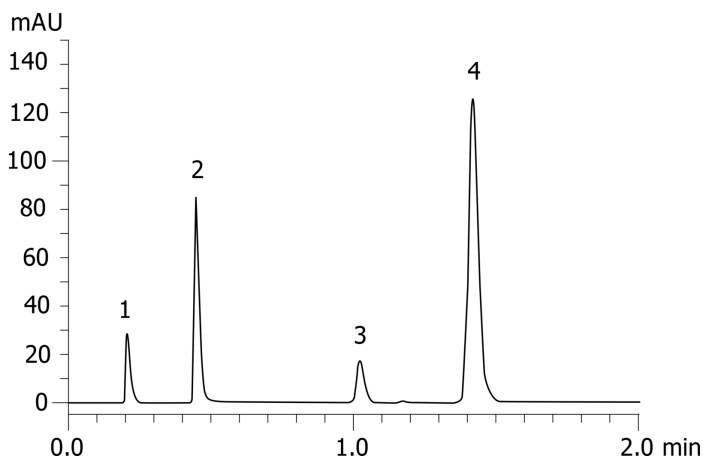

### ANALYTES:

- 1** Uracil  
Retention Time: 0.199 min
- 2** Acetophenone  
Retention Time: 0.435 min
- 3** Toluene  
Retention Time: 0.994 min
- 4** Naphthalene  
Retention Time: 1.38 min

