

# HPLC Application

ID No.: 19980

## Kinetex 2.6 µm XB-C18 50 x 4.6 Prodigy Test With SecurityGuard ULTRA

**Column:** Kinetex® 2.6 µm XB-C18 100 Å, LC Column 50 x 4.6 mm, Ea

**Dimensions:** 50 x 4.6 mm ID

**Order No:** 00B-4496-E0

**Elution Type:** Isocratic

**Eluent A:** Water

**Eluent B:** Acetonitrile

Gradient Profile:	Step No.	Time (min)	Pct A	Pct B
	1	0	35	65


**Flow Rate:** 1.85 mL/min

**Col. Temp.:** ambient

**Detection:** UV-Vis Abs.-Diode Array (PDA) @ (ambient)  
19980



Products used in this application:

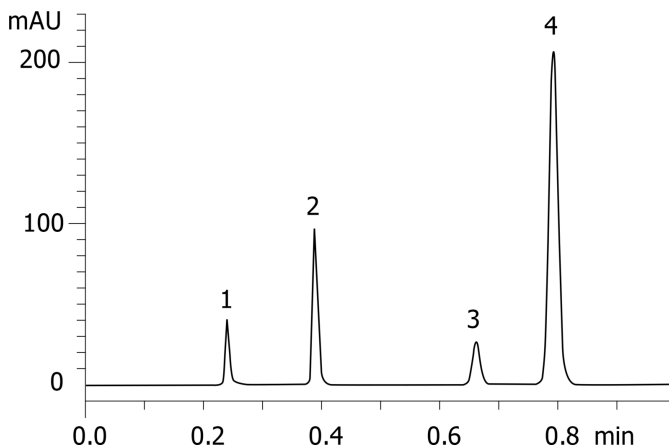

### ANALYTES:

- 1** Uracil  
Retention Time: 0.24 min
- 2** Acetophenone  
Retention Time: 0.39 min
- 3** Toluene  
Retention Time: 0.661 min
- 4** Naphthalene  
Retention Time: 0.793 min

