

## Overlay: Yarra 3µm SEC-3000 and BioSep-SEC-S 3000 Low Salt PMQC Standard

**Column:** Yarra™ 3 µm SEC-3000, LC Column 300 x 7.8 mm, Ea

**Dimensions:** 300 x 7.8 mm ID

**Order No:** 00H-4513-K0

**Elution Type:** Isocratic


**Eluent A:** 100 mM Sodium Phosphate in Water pH 6.8 with 0.025% Na3

Gradient Profile:	Step No.	Time (min)	Pct A
	1	0	100

**Flow Rate:** 1 mL/min

**Col. Temp.:** ambient

**Detection:** UV-Vis Abs.-Variable Wave.(UV) @ 280 nm (ambient)

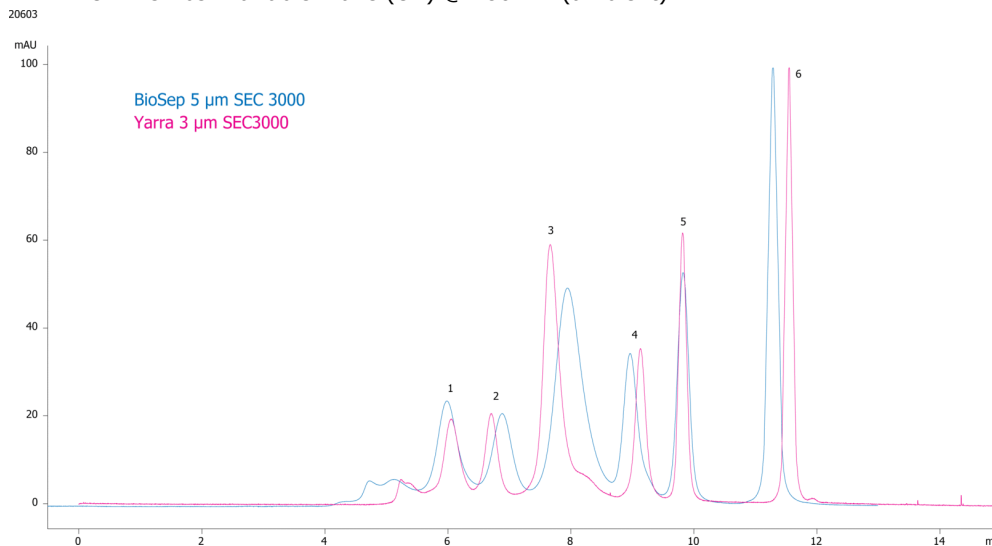

Products used in this application:



### ANALYTES:

- 1** Thyroglobulin  
Retention Time: 6 min
- 2** IgA  
Retention Time: 6.7 min
- 3** IgG  
Retention Time: 7.7 min
- 4** Ovalbumin  
Retention Time: 8.2 min
- 5** Myoglobin  
Retention Time: 9.8 min
- 6** Uridine  
Retention Time: 11 min

