

## Yarra 3µm SEC-2000 low salt buffer PMQC Standard

**Column:** Yarra™ 3 µm SEC-2000, LC Column 300 x 7.8 mm, Ea

**Dimensions:** 300 x 7.8 mm ID

**Order No:** 00H-4512-K0

**Elution Type:** Isocratic

**Eluent A:** 100 mM Sodium Phosphate in Water pH 6.8 with 0.025% Na3


Gradient	Step No.	Time (min)	Pct A
<b>Profile:</b>	<b>1</b>	<b>0</b>	<b>100</b>

**Flow Rate:** 1 mL/min

**Col. Temp.:** ambient

**Detection:** UV-Vis Abs.-Variable Wave.(UV) @ 280 nm (ambient)

20613

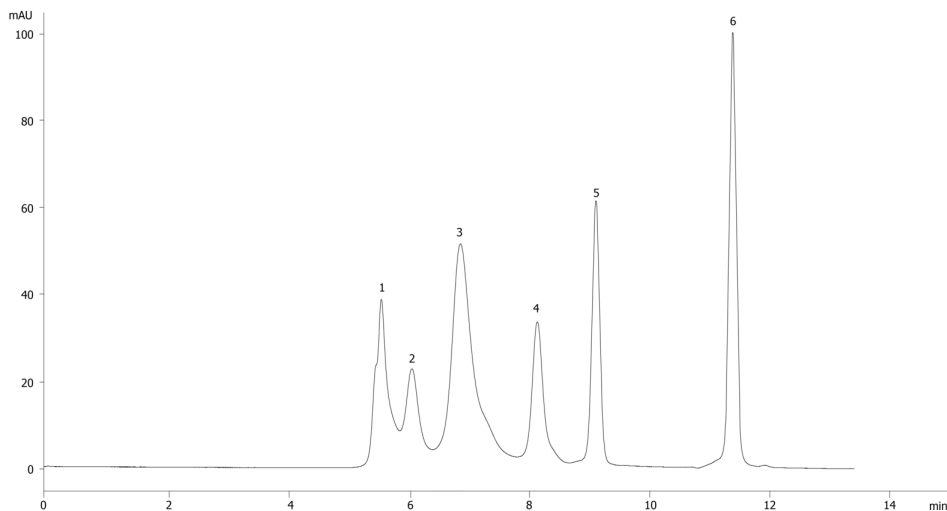

### ANALYTES:

- 1** Thyroglobulin  
Retention Time: 5.4 min
- 2** IgA  
Retention Time: 5.9 min
- 3** IgG  
Retention Time: 6.7 min
- 4** Ovalbumin  
Retention Time: 8.2 min
- 5** Myoglobin  
Retention Time: 9 min
- 6** Uridine  
Retention Time: 11.3 min



Products used in this application:

