

Overlay: Yarra 3µm SEC-2000 vs TSKgel super SW2000 low salt buffer PMQC Standard

Column: Yarra™ 3 µm SEC-2000, LC Column 300 x 7.8 mm, Ea

Dimensions: 300 x 7.8 mm ID

Order No: 00H-4512-K0

Elution Type: Isocratic

Eluent A: 100 mM Sodium Phosphate in Water pH 6.8 with 0.025% Na3


Gradient Profile:	Step No.	Time (min)	Pct A
	1	0	100

Flow Rate: 1 mL/min

Col. Temp.: ambient

Detection: UV-Vis Abs.-Variable Wave.(UV) @ 280 nm (ambient)

20615

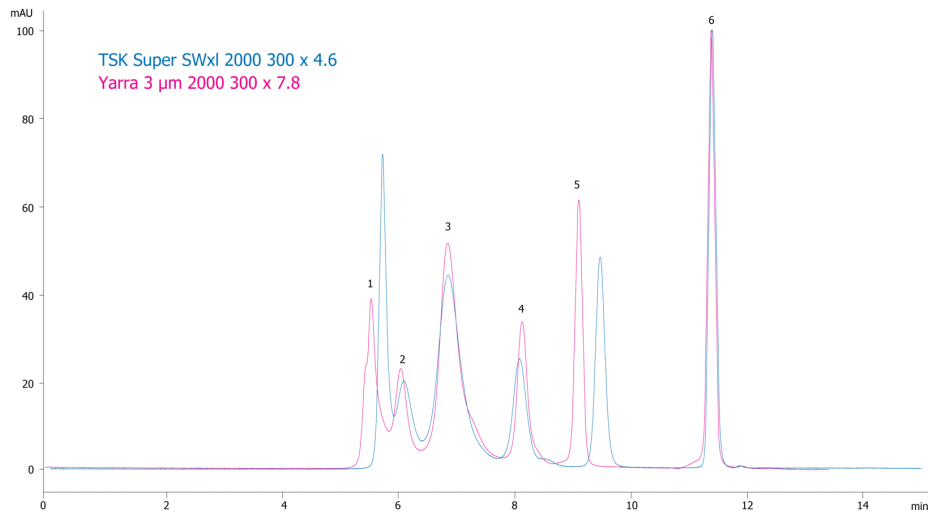

Products used in this application:



ANALYTES:

- 1** Thyroglobulin
Retention Time: 5.5 min
- 2** IgA
Retention Time: 6.2 min
- 3** IgG
Retention Time: 7 min
- 4** Ovalbumin
Retention Time: 8.2 min
- 5** Myoglobin
Retention Time: 9 min
- 6** Uridine
Retention Time: 11.5 min

