

## Calibration curves for Yarra phases SEC-2000, SEC-3000, SEC-4000

**Column:** Yarra™ 3 µm SEC-2000, LC Column 300 x 7.8 mm, Ea

**Dimensions:** 300 x 7.8 mm ID

**Order No:** 00H-4512-K0

**Elution Type:** Isocratic

**Eluent A:** 100 mM Sodium Chloride pH 6.8

Gradient Profile:	Step No.	Time (min)	Pct A
	1	0	100

**Flow Rate:** 1 mL/min


**Col. Temp.:** ambient

**Detection:** UV-Vis Abs.-Variable Wave.(UV) @ 280 nm (ambient)

20638



Products used in this application:

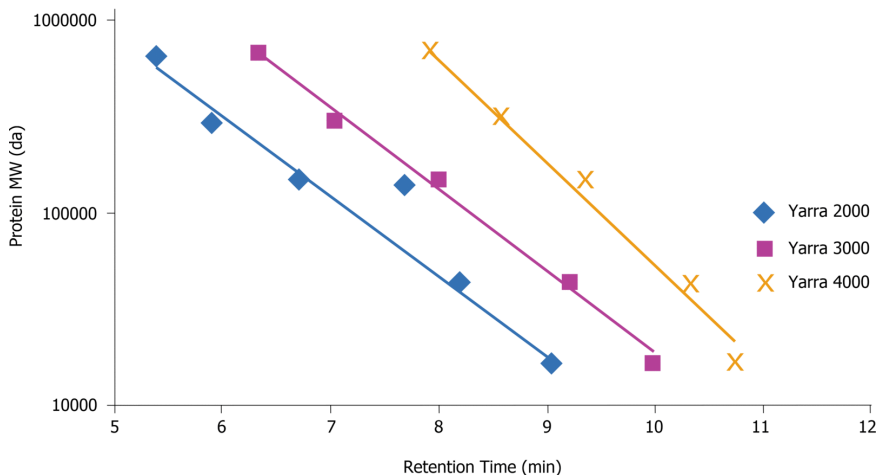

### ANALYTES:

- 1 Thyroglobulin
- 2 IgA
- 3 IgG
- 4 Ovalbumin
- 5 Myoglobin
- 6 Uridine

