

## Degraded Myoglobin on Aeris WIDEPORE 3.6u XB-C18 vs. Jup 5u C18 overlay

**Column:** Aeris™ 3.6 µm WIDEPORE XB-C18 200 Å, LC Column 150 x 2.1 mm, Ea

**Dimensions:** 150 x 2.1 mm ID

**Order No:** 00F-4482-AN

**Elution Type:** Gradient

**Eluent A:** 0.1% TFA

**Eluent B:** Acetonitrile

Gradient Profile:	Step No.	Time (min)	Pct A	Pct B
	1	0	90	10
	2	5	90	10
	3	50	50	50

**Flow Rate:** 0.2 mL/min


**Col. Temp.:** ambient

**Detection:** UV-Vis Abs.-Diode Array (PDA) @ (ambient)

**Analyst Note:** SecurityGuard™ Guard Cartridge System extends column lifetime.

- SecurityGuard ULTRA Cartridges, UHPLC WIDEPORE C18 for 2.1mm ID, 3/Pk Part No.: AJ0-8783


- Holder Part No.: AJ0-9000

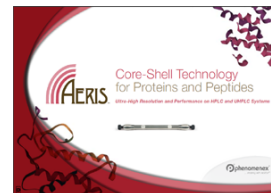

Products used in this application:



20682

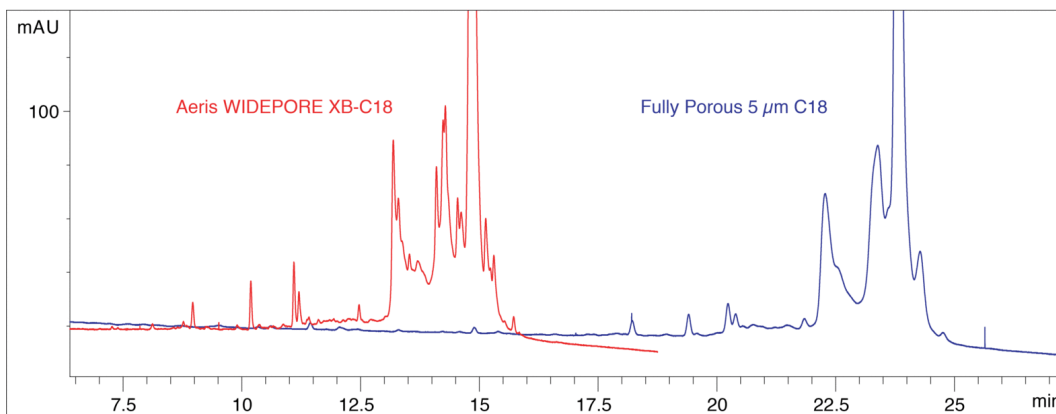

### ANALYTES:

1 Myoglobin

