

## Yarra 3µm SEC-3000 AcCN gradient wash Ig-G + BSA before and after overlay

**Column:** Yarra™ 3 µm SEC-3000, LC Column 300 x 7.8 mm, Ea

**Dimensions:** 300 x 7.8 mm ID

**Order No:** 00H-4513-K0

**Elution Type:** Isocratic

**Eluent A:** 100 mM NaPO4 pH 6.8



Gradient Profile:	Step No.	Time (min)	Pct A
	1	0	100

**Flow Rate:** 1 mL/min

**Col. Temp.:** ambient

**Detection:** UV-Vis Abs.-Variable Wave.(UV) @ 280 nm (ambient)

20744

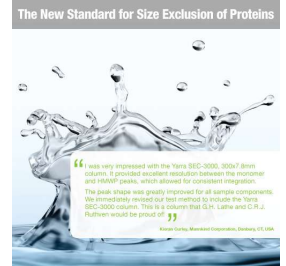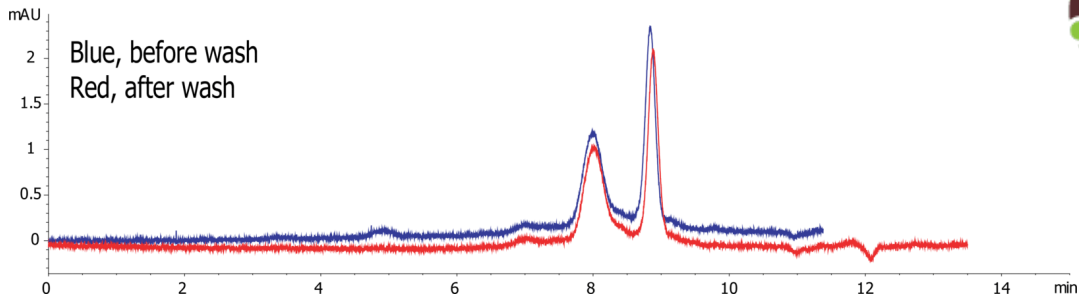

Products used in this application:



### ANALYTES:

- 1 IgG monomer  
Retention Time: 8 min
- 2 Albumin  
Retention Time: 8.9 min