

HPLC Application

ID No.: 20901

Quantification of Acrylamide in Food by LC/MS/MS

Column: Luna® 3 µm C18(2) 100 Å, LC Column 150 x 3 mm, Ea

Dimensions: 150 x 3 mm ID

Order No: 00F-4251-Y0

Elution Type: Gradient

Eluent A: Water with 0.001% Formic Acid

Eluent B: Methanol with 0.001% Formic Acid

Gradient Profile:	Step No.	Time (min)	Pct A	Pct B
	1	0	90	10
	2	10	100	0
	3	0	90	10

Flow Rate: 250 µL/min


Col. Temp.: 60 °C

Detection: Tandem Mass Spec (MS-MS) @ (ambient)

Detector Info: SCIEX
20901



Products used in this application:

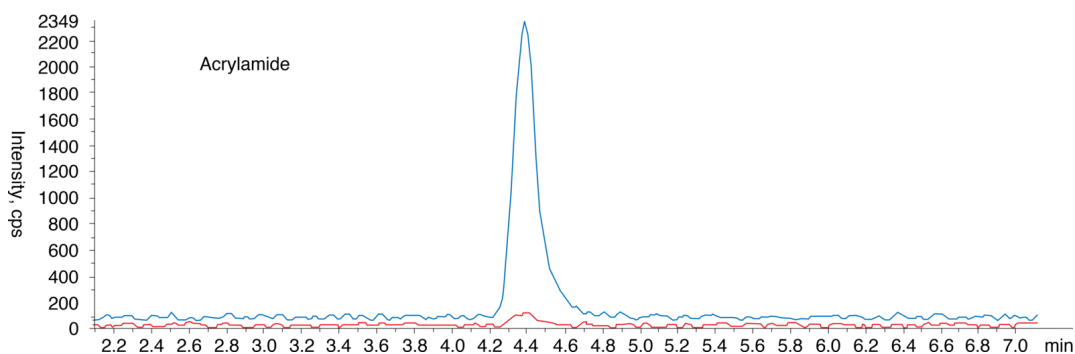

ANALYTES:

- 1 Acrylamide
Retention Time: 4.4 min

