

## Overlay: Yarra 3µm SEC-2000 with different NaSO<sub>4</sub> MP


**Column:** Yarra™ 3 µm SEC-2000, LC Column 300 x 7.8 mm, Ea  
**Dimensions:** 300 x 7.8 mm ID  
**Order No:** 00H-4512-K0  
**Elution Type:** Isocratic  
**Eluent A:** 50 mM Sodium Phosphate in Water pH 6.8 + various salt concentrations  
**Gradient**  

Step No.	Time (min)	Pct A
1	0	100

  
**Flow Rate:** 1 mL/min  
**Col. Temp.:** ambient  
**Detection:** UV-Vis Abs.-Variable Wave.(UV) @ 280 nm (ambient)  
21464



Products used in this application:

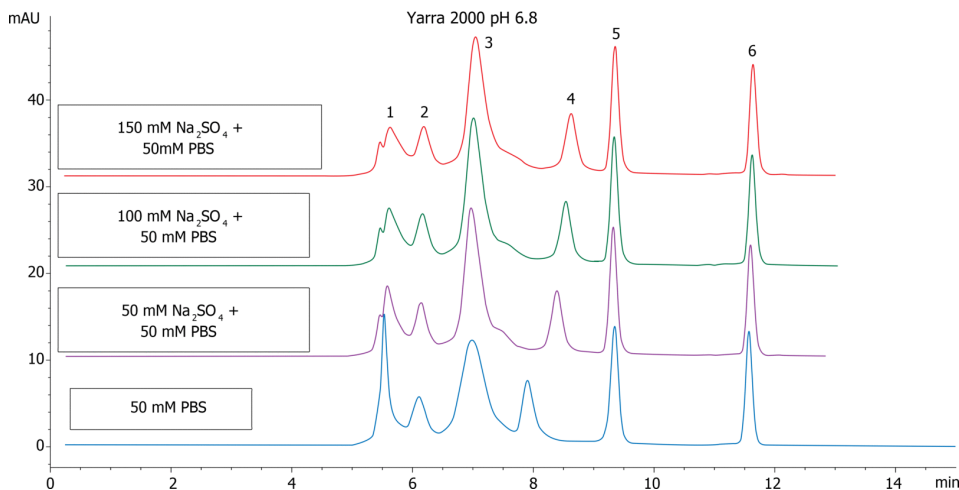

### ANALYTES:

- 1** Thyroglobulin  
Retention Time: 5.4 min
- 2** IgA  
Retention Time: 5.9 min
- 3** IgG  
Retention Time: 6.7 min
- 4** Ovalbumin  
Retention Time: 8.2 min
- 5** Myoglobin  
Retention Time: 9 min
- 6** Uridine  
Retention Time: 11.3 min

