

# HPLC Application

ID No.: 21570

## FMOC-Tyr(tBu)-OH on Lux 5 $\mu$ m Cellulose-3 in RP

**Column:** Lux<sup>®</sup> 5  $\mu$ m Cellulose-3, LC Column 250 x 4.6 mm, Ea

**Dimensions:** 250 x 4.6 mm ID

**Order No:** 00G-4493-E0

**Elution Type:** Isocratic

**Eluent A:** ACN / 0.1% TFA (60:40) TFA (0.1%)


Gradient	Step No.	Time (min)	Pct A
<b>Profile:</b>	<b>1</b>	0	100

**Flow Rate:** 1 mL/min

**Col. Temp.:** ambient

**Detection:** UV-Vis Abs.-Variable Wave.(UV) @ 220 nm (ambient)

21570

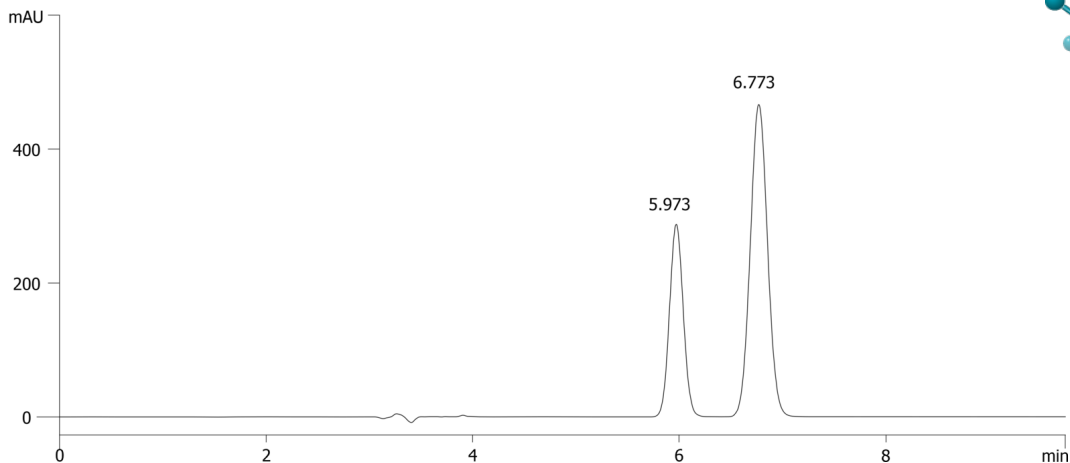

Products used in this application:



### ANALYTES:

- 1** FMOC-Tyr(tBu)-OH  
Retention Time: 5.973 min
- 2** FMOC-Tyr(tBu)-OH  
Retention Time: 6.773 min

