

## FMOC-His(Trt)-OH on Lux 5µm Cellulose-1 in RP

**Column:** Lux® 5 µm Cellulose-1, LC Column 250 x 4.6 mm, Ea

**Dimensions:** 250 x 4.6 mm ID

**Order No:** 00G-4459-E0

**Elution Type:** Isocratic

**Eluent A:** MeOH / 0.1% Formic Acid (90:10) Formic acid (0.1%)


| Gradient        | Step No. | Time (min) | Pct A |
|-----------------|----------|------------|-------|
| <b>Profile:</b> | <b>1</b> | 0          | 100   |

**Flow Rate:** 1 mL/min

**Col. Temp.:** ambient

**Detection:** UV-Vis Abs.-Variable Wave.(UV) @ 220 nm (ambient)

21583

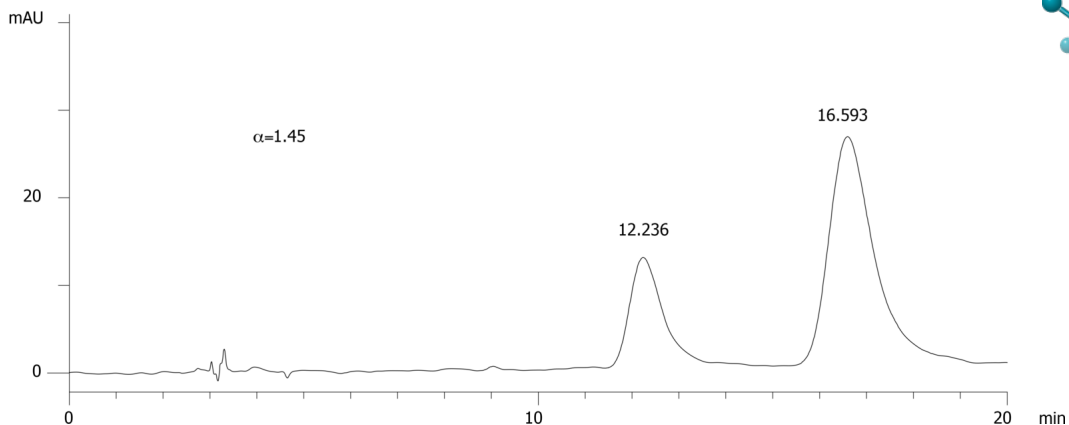

Products used in this application:



### ANALYTES:

- 1** FMOC-His(Trt)-OH  
Retention Time: 12.236 min
- 2** FMOC-His(Trt)-OH  
Retention Time: 16.593 min

