

## FMOC-Val-OH on Lux 5µm Cellulose-3 in RP

**Column:** Lux® 5 µm Cellulose-3, LC Column 250 x 4.6 mm, Ea

**Dimensions:** 250 x 4.6 mm ID

**Order No:** 00G-4493-E0

**Elution Type:** Isocratic

**Eluent A:** MeOH / 0.1% Formic Acid (80:20) with None in buffer

Gradient	Step No.	Time (min)	Pct A
Profile:	1	0	100


**Flow Rate:** 1 mL/min

**Col. Temp.:** ambient

**Detection:** UV-Vis Abs.-Variable Wave.(UV) @ 220 nm (ambient)  
21657



Products used in this application:

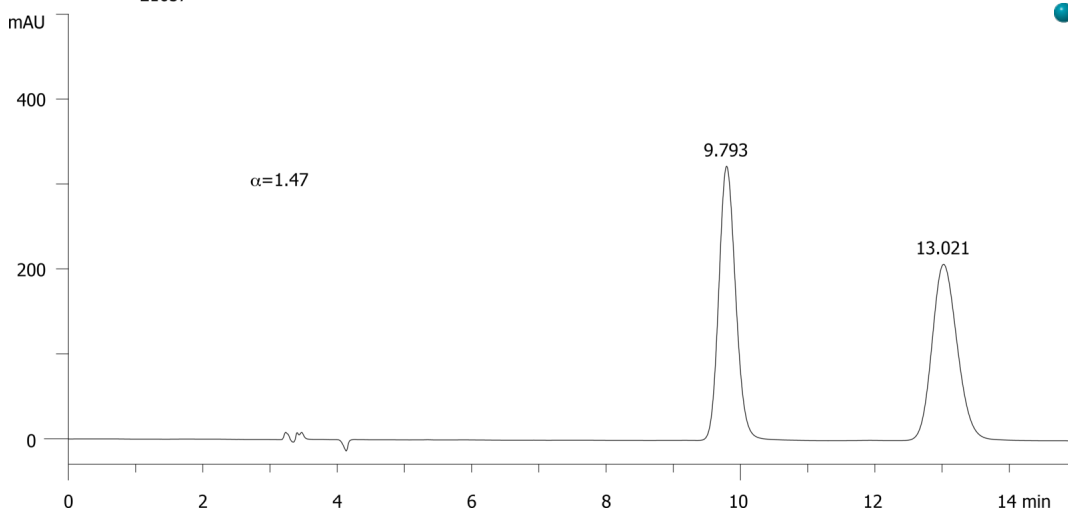

### ANALYTES:

1 FMOC-Val-OH

Retention Time: 9.793 min

2 FMOC-Val-OH

Retention Time: 13.021 min

