

## Cyproconazole on Lux 5µm Cellulose-2 in RP


**Column:** Lux® 5 µm Cellulose-2, LC Column 250 x 4.6 mm, Ea  
**Dimensions:** 250 x 4.6 mm ID  
**Order No:** 00G-4457-E0  
**Elution Type:** Isocratic  
**Eluent A:** ACN / 20mM NH<sub>4</sub>HCO<sub>3</sub>+0.1%DEA (60:40) with DEA (0.1%) in buffer  
**Gradient**

Step No.	Time (min)	Pct A
1	0	100

**Profile:**  
**Flow Rate:** 1 mL/min  
**Col. Temp.:** ambient  
**Detection:** UV-Vis Abs.-Variable Wave.(UV) @ 220 nm (ambient)  
App ID 21754



Products used in this application:

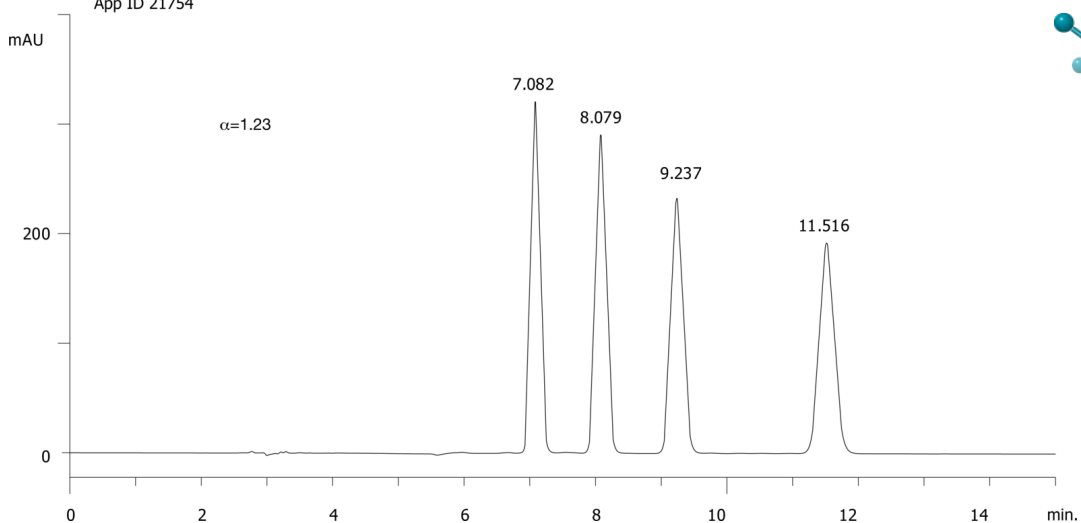

### ANALYTES:

- 1 Cyproconazole  
Retention Time: 7.082 min
- 2 Cyproconazole  
Retention Time: 8.079 min
- 3 Cyproconazole  
Retention Time: 9.237 min
- 4 Cyproconazole  
Retention Time: 11.516 min

