

## Cyproconazole on Lux 5 $\mu$ m Cellulose-2 in RP


**Column:** Lux<sup>®</sup> 5  $\mu$ m Cellulose-2, LC Column 250 x 4.6 mm, Ea  
**Dimensions:** 250 x 4.6 mm ID  
**Order No:** 00G-4457-E0  
**Elution Type:** Isocratic  
**Eluent A:** MeOH / 20mM NH<sub>4</sub>HCO<sub>3</sub>+0.1%DEA (80:20) with DEA (0.1%) in buffer  
**Gradient**

Step No.	Time (min)	Pct A
1	0	100

**Profile:**  
**Flow Rate:** 1 mL/min  
**Col. Temp.:** ambient  
**Detection:** UV-Vis Abs.-Variable Wave.(UV) @ 220 nm (ambient)  
App ID 21756



Products used in this application:

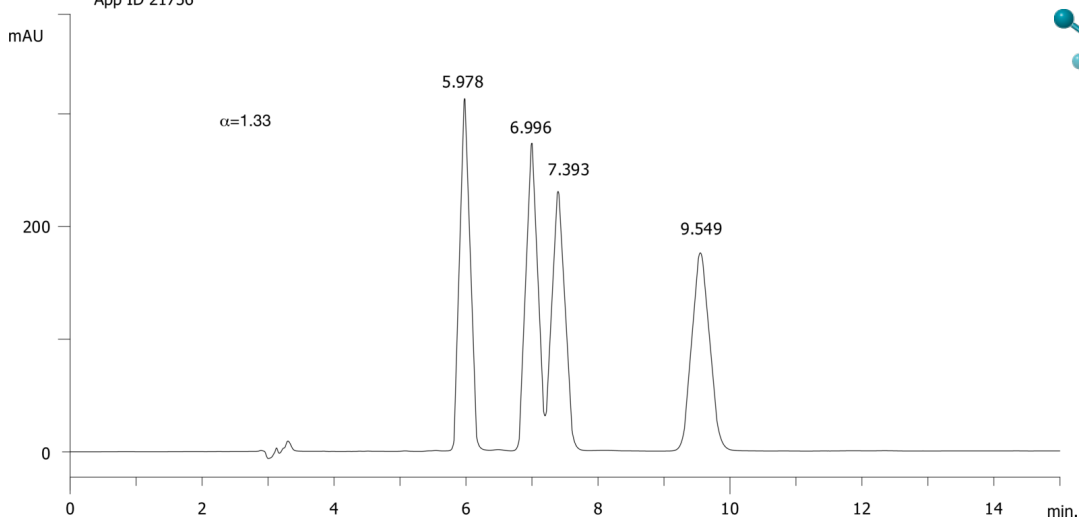

### ANALYTES:

- 1 Cyproconazole  
Retention Time: 5.978 min
- 2 Cyproconazole  
Retention Time: 6.996 min
- 3 Cyproconazole  
Retention Time: 7.393 min
- 4 Cyproconazole  
Retention Time: 9.549 min

