

# HPLC Application

ID No.: 22121

## Fast Pain Panel using Kinetex 2.6 u Biphenyl 50x3.0mm

**Column:** Kinetex<sup>®</sup> 2.6  $\mu$ m Biphenyl 100  $\text{\AA}$ , LC Column 50 x 3.0 mm, Ea

**Dimensions:** 50 x 3 mm ID

**Order No:** 00B-4622-Y0

**Elution Type:** Gradient

**Eluent A:** Water with 0.1% formic acid

**Eluent B:** Methanol with 0.1% formic acid

Gradient Profile:	Step No.	Time (min)	Pct A	Pct B
	1	0	90	10
	2	2.5	0	100
	3	3.5	0	100
	4	3.51	90	10
	5	5	90	10

**Flow Rate:** 0.7 mL/min


**Col. Temp.:** 40 °C

**Detection:** Tandem Mass Spec (MS-MS) @ (40 °C)

22121



Products used in this application:

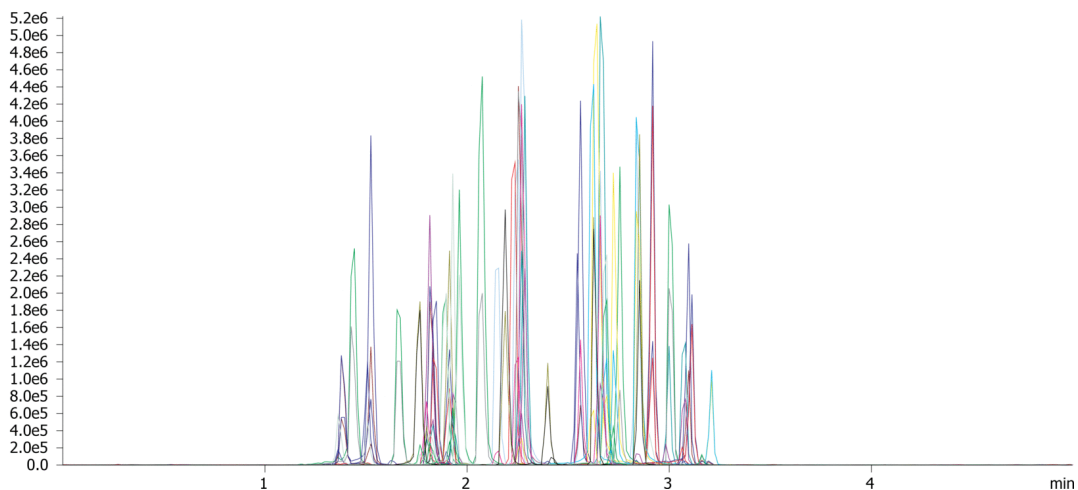

## Fast Pain Panel using Kinetex 2.6 u Biphenyl 50x3.0mm

### ANALYTES:

- |                             |                    |
|-----------------------------|--------------------|
| <b>1</b> Morphine           | <b>41</b> Diazepam |
| <b>2</b> Oxymorphone        |                    |
| <b>3</b> Hydromorphone      |                    |
| <b>4</b> Amphetamine        |                    |
| <b>5</b> Naloxone           |                    |
| <b>6</b> Methamphetamine    |                    |
| <b>7</b> Codeine            |                    |
| <b>8</b> MDA                |                    |
| <b>9</b> Oxycodone          |                    |
| <b>10</b> Naltrexone        |                    |
| <b>11</b> Hydrocodone       |                    |
| <b>12</b> MDMA              |                    |
| <b>13</b> MDEA              |                    |
| <b>14</b> Norfentanyl       |                    |
| <b>15</b> Tramadol          |                    |
| <b>16</b> Benzoylcegonine   |                    |
| <b>17</b> Meperidine        |                    |
| <b>18</b> Meprobamate       |                    |
| <b>19</b> Norbuprenorphine  |                    |
| <b>20</b> Fentanyl          |                    |
| <b>21</b> Buprenorphine     |                    |
| <b>22</b> Flurazepam        |                    |
| <b>23</b> Carisoprodol      |                    |
| <b>24</b> PCP               |                    |
| <b>25</b> Propoxyphene      |                    |
| <b>26</b> Sufentanil        |                    |
| <b>27</b> 6-MAM             |                    |
| <b>28</b> Midazolam         |                    |
| <b>29</b> Normeperidine     |                    |
| <b>30</b> EDDP              |                    |
| <b>31</b> Methadone         |                    |
| <b>32</b> Lorazepam         |                    |
| <b>33</b> Clonazepam        |                    |
| <b>34</b> Norpropoxyphene   |                    |
| <b>35</b> Oxazepam          |                    |
| <b>36</b> Hydroxyalprazolam |                    |
| <b>37</b> Nordiazepam       |                    |
| <b>38</b> Flunitrazepam     |                    |
| <b>39</b> Temazepam         |                    |
| <b>40</b> Alprazolam        |                    |

