

## Apo-Myoglobin Tryptic Digest on Aeris PEPTIDE 3.6u XB-C18 150x2.1mm

**Column:** Aeris™ 3.6 µm PEPTIDE XB-C18 100 Å, LC Column 150 x 2.1 mm, Ea

**Dimensions:** 150 x 2.1 mm ID

**Order No:** 00F-4507-AN

**Elution Type:** Gradient

**Eluent A:** 0.1% Formic Acid in Water

**Eluent B:** 0.1% Formic Acid in Acetonitrile


| Gradient Profile: | Step No. | Time (min) | Pct A | Pct B |
|-------------------|----------|------------|-------|-------|
|                   | 1        | 0          | 97    | 3     |
|                   | 2        | 1.5        | 97    | 3     |
|                   | 3        | 15         | 35    | 65    |
|                   | 4        | 17         | 5     | 95    |

**Flow Rate:** 0.5 mL/min


**Col. Temp.:** ambient

**Detection:** UV-Vis Abs.-Variable Wave.(UV) @ 214 nm (ambient)

22136

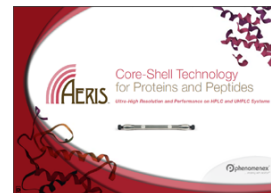

Products used in this application:

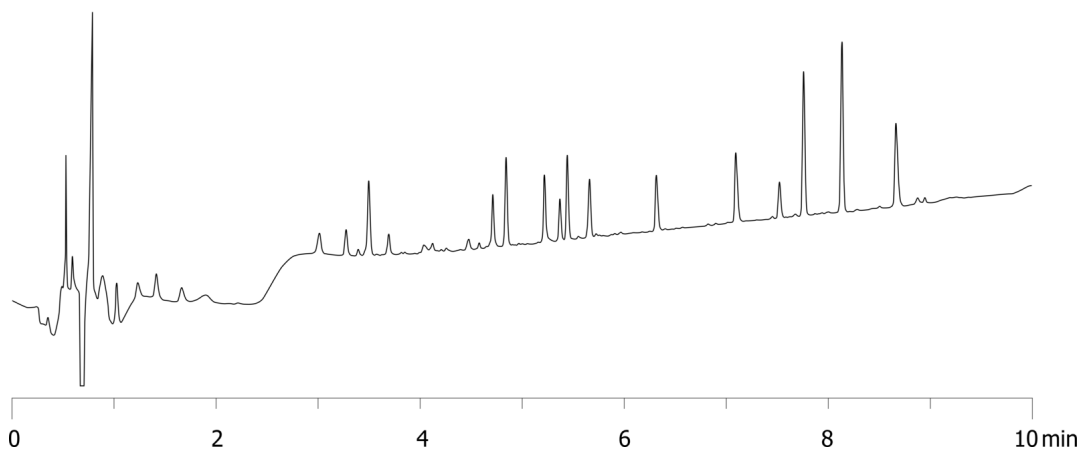

### ANALYTES:

- 1 Apo Myoglobin Tryptic Digest

