

# HPLC Application

ID No.: 22137

## Apo-Myoglobin Tryptic Digest on Aeris PEPTIDE 1.7u XB-C18 150x2.1mm

**Column:** Aeris<sup>™</sup> 1.7  $\mu$ m PEPTIDE XB-C18 100  $\text{\AA}$ , LC Column 150 x 2.1 mm, Ea

**Dimensions:** 150 x 2.1 mm ID

**Order No:** 00F-4506-AN

**Elution Type:** Gradient

**Eluent A:** 0.1% Formic Acid in Water

**Eluent B:** 0.1% Formic Acid in Acetonitrile


Gradient Profile:	Step No.	Time (min)	Pct A	Pct B
	1	0	97	3
	2	1.5	97	3
	3	15	35	65
	4	17	5	95

**Flow Rate:** 0.5 mL/min


**Col. Temp.:** ambient

**Detection:** UV-Vis Abs.-Variable Wave.(UV) @ 214 nm (ambient)

22137

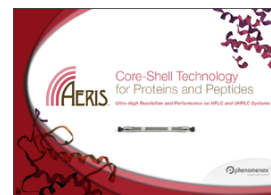

Products used in this application:

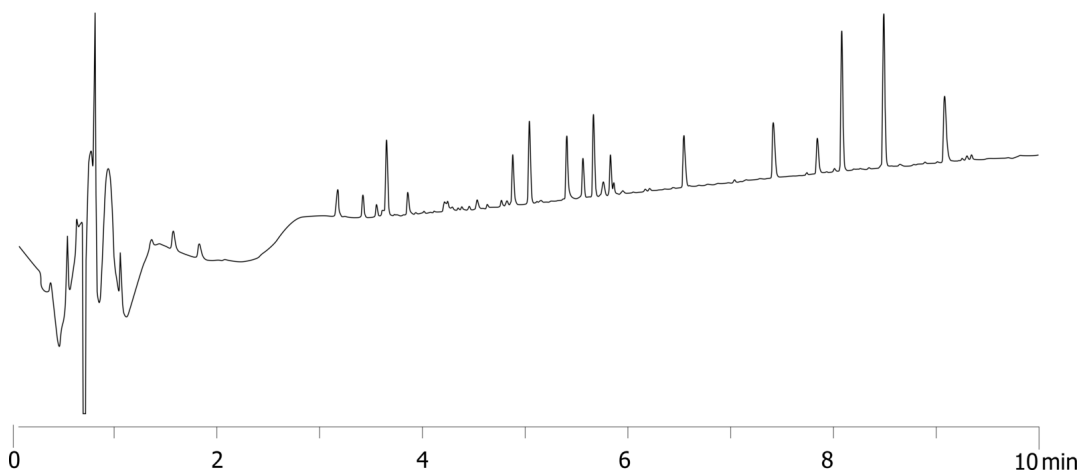

### ANALYTES:

- 1 Apo Myoglobin Tryptic Digest

