

HPLC Application

ID No.: 22140

Apo-Myoglobin Tryptic Digest on Aeris PEPTIDE 3.6u XB-C18 250x2.1mm FAST GRADIENT

Column: Aeris[™] 3.6 μ m PEPTIDE XB-C18 100 \AA , LC Column 250 x 4.6 mm, Ea

Dimensions: 250 x 4.6 mm ID

Order No: 00G-4507-E0

Elution Type: Gradient

Eluent A: 0.1% Formic Acid in Water

Eluent B: 0.1% Formic Acid in Acetonitrile


Gradient Profile:	Step No.	Time (min)	Pct A	Pct B
	1	0	97	3
	2	1	97	3
	3	10	35	65
	4	12.5	5	95

Flow Rate: 0.9 mL/min


Col. Temp.: ambient

Detection: UV-Vis Abs.-Variable Wave.(UV) @ 214 nm (ambient)

22140

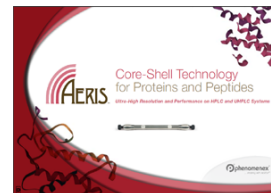

Products used in this application:

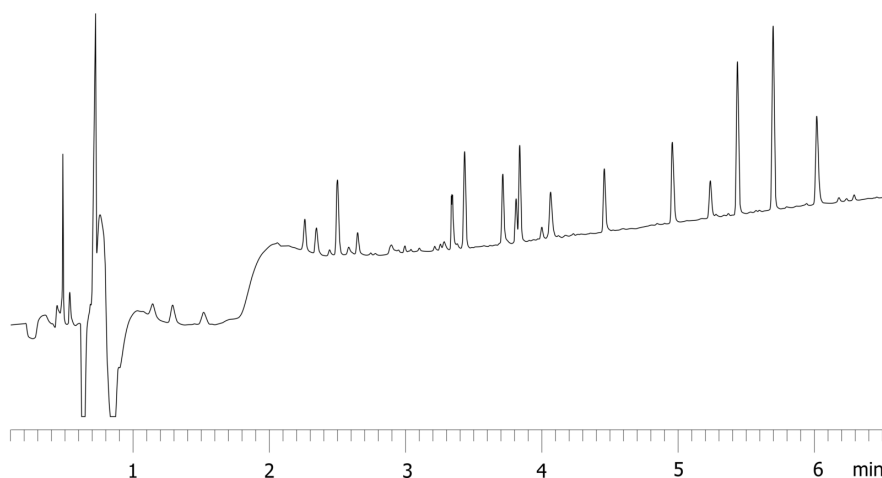

ANALYTES:

- 1 Apo Myoglobin Tryptic Digest

