

# HPLC Application

ID No.: 22145

## Steroid Panel on Kinetex 2.6u Biphenyl 50x3.0

**Column:** Kinetex<sup>®</sup> 2.6  $\mu$ m Biphenyl 100  $\text{\AA}$ , LC Column 50 x 3.0 mm, Ea

**Dimensions:** 50 x 3 mm ID

**Order No:** 00B-4622-Y0

**Elution Type:** Gradient

**Eluent A:** Water

**Eluent B:** Acetonitrile

Gradient Profile:	Step No.	Time (min)	Pct A	Pct B
	1	0	80	20
	2	5	40	60
	3	6	40	60


**Flow Rate:** 0.6 mL/min

**Col. Temp.:** 25  $^{\circ}$ C

**Detection:** UV-Vis Abs.-Variable Wave.(UV) @ (ambient)  
22145



Products used in this application:

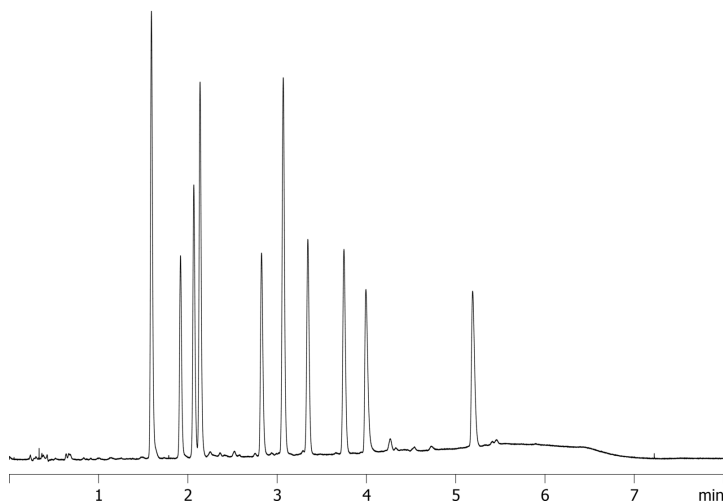

## Steroid Panel on Kinetex 2.6u Biphenyl 50x3.0

### ANALYTES:

- 1 Estriol
- 2 Hydrocortisone
- 3 Predisone
- 4 Cortisone
- 5 Corticosterone
- 6 Beta-Estradiol
- 7 Cortisone Acetate
- 8 21-Hydroxyprogesterone
- 9 17-Hydroxyprogesterone
- 10 Deoxycorticosterone

