# **HPLC Application** ID No.: **22770**



### Oxytocin and ARG-Vasopressin Chromatogram

Luna® 3 µm PFP(2) 100 Å, LC Column 50 x 2 mm, Ea

50 x 2 mm ID **Dimensions:** Order No: 00B-4447-B0 Elution Type: Gradient

Eluent A: 0.1%FA in water

Eluent B: 0.1%FA in 50:50 MeOH:Acetonitrile

Gradient	Step No.	Time (min)	Pct A	Pct B
Profile:	1	0	90	10
	2	0.2	90	10
	3	0.5	70	30
	4	3.5	70	30
	5	3.51	30	70
	6	4	30	70
	7	4.01	90	10
	8	5	90	10




Products used in this application:



Flow Rate: 0.45 mL/min Col. Temp.: ambient

**Detection:** Tandem Mass Spec (MS-MS) @ (40 °C)

**Detector Info:** <a target="\_blank"

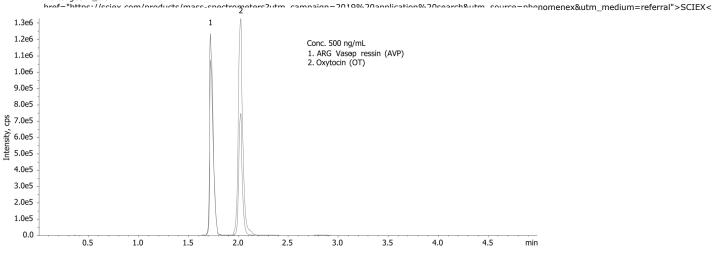

### **ANALYTES:**

1 Oxytocin

Retention Time: 2.02 min

ARG-Vasopressin

Retention Time: 1.72 min

©2025 Phenomenex Inc. All rights reserved.

For more information contact your Phenomenex Representative at support@phenomenex.com



Phenomenex products are available worldwide. www.phenomenex.com

## Sample Preparation Details for HPLC Application ID No.: 22770



## **Oxytocin and ARG-Vasopressin Chromatogram**

#### **PRODUCT DESCRIPTION:**

Strata™-X-CW 33 µm Polymeric Weak Cation, 30 mg / well, 96-Well Plates , 2/Pk

Order No.: 8E-S035-TGB

SOLID PHASE EXTRACTION (SPI	E) PRODCEDURE
-----------------------------	---------------

Note: The solvent volumes shown below are for a 30 mg bed mass.

The solvent volumes will need to be adjusted for a smaller or larger bed mass.

Condition:	
Load:	
Wash:	
Dry:	
3-5 min	
Elute:	
Final Prep and Analysis:	
Inject: 10 μL on HPLC Tandem Mass Spec (MS-MS) @ (40°C)	

ANALYTES:	Spiked Conc.	Log P	рКа	% Rec	%RSC
	(ng/mL)				(n=0)
<b>1</b> Oxytocin	0				
2 ARG-Vasopressin	0				

This method is designed as a convenient starting point for further investigation and can be tailored to meet your extraction goals. Call your local Phenomenex Representative for assistance in method development and optimization techniques.

©2025 Phenomenex Inc. All rights reserved.

For more information contact your Phenomenex Representative at support@phenomenex.com

