

HPLC Application

ID No.: 23055

Fast Analysis of Sugars in Fruit Juices - Merlot (Red Wine)

Column: Rezex[™] RCM-Monosaccharide Ca+2 (8%), LC Column 100 x 7.8 mm, Ea

Dimensions: 100 x 7.8 mm ID

Order No: 00D-0130-K0

Elution Type: Isocratic

Eluent A: Water

| Gradient | Step No. | Time (min) | Pct A |
|-----------------|----------|------------|-------|
| Profile: | 1 | 0 | 100 |

Flow Rate: 0.4 mL/min

Col. Temp.: 80 °C


Detection: Refractive Index (RI) @ (35 °C)
23055

REZEX[™]
HPLC COLUMNS for Carbohydrates, Oligosaccharides, and Organic Acid Analysis



Products used in this application:

REZEX[™]

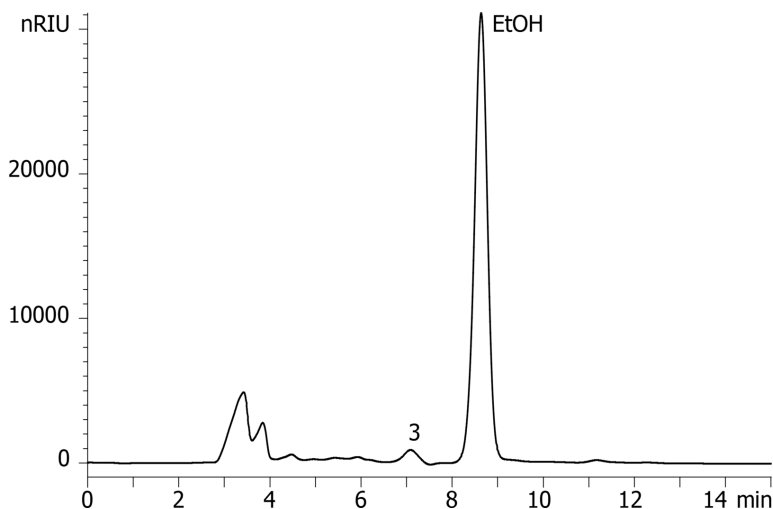

ANALYTES:

1 Merlot (Red Wine)

