

# HPLC Application

ID No.: 23059

## Fast Analysis of Sugars in Fruit Juices - Fruit Punch

**Column:** Rezex<sup>™</sup> RCM-Monosaccharide Ca+2 (8%), LC Column 100 x 7.8 mm, Ea

**Dimensions:** 100 x 7.8 mm ID

**Order No:** 00D-0130-K0

**Elution Type:** Isocratic

**Eluent A:** Water

Gradient	Step No.	Time (min)	Pct A
<b>Profile:</b>	<b>1</b>	<b>0</b>	<b>100</b>

**Flow Rate:** 0.4 mL/min

**Col. Temp.:** 80 °C


**Detection:** Refractive Index (RI) @ (35 °C)  
23059

**REZEX**<sup>™</sup>  
HPLC COLUMNS for Carbohydrates, Oligosaccharides, and Organic Acid Analysis



Products used in this application:

**REZEX**<sup>™</sup>

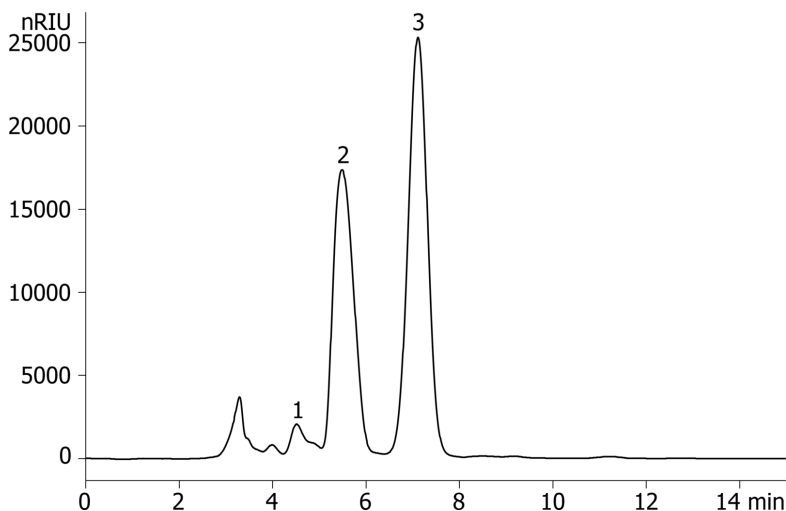

### ANALYTES:

1 Fruit Punch

