

HPLC Application

ID No.: 24416

BEH Comps - Pesticide Panel using Luna Omega 1.6um Polar C18 50 x 2.1mm

Column: Luna Omega 1.6 um Polar C18 100 A LC Column 50 x 2.1 mm, Ea

Dimensions: 50 x 2.1 mm ID

Order No: 00B-4748-AN

Elution Type: Gradient

Eluent A: 10mM Ammonium formate

Eluent B: 0.1% Formic acid in Acetonitrile

Gradient Profile:	Step No.	Time (min)	Pct A	Pct B
	1	0	97	3
	2	5	0	100
	3	7	0	100

Flow Rate: 0.5 mL/min

Col. Temp.: 40 °C

Detection: Mass Spectrometer (MS) @ amu (ambient)


Detector Info: <a target="_blank"

24416

href="https://sciex.com/products/mass-spectrometers?utm_campaign=2019%20application%20search&utm_source=phenomenex&utm_medium=referral">SCIEX<



Products used in this application:

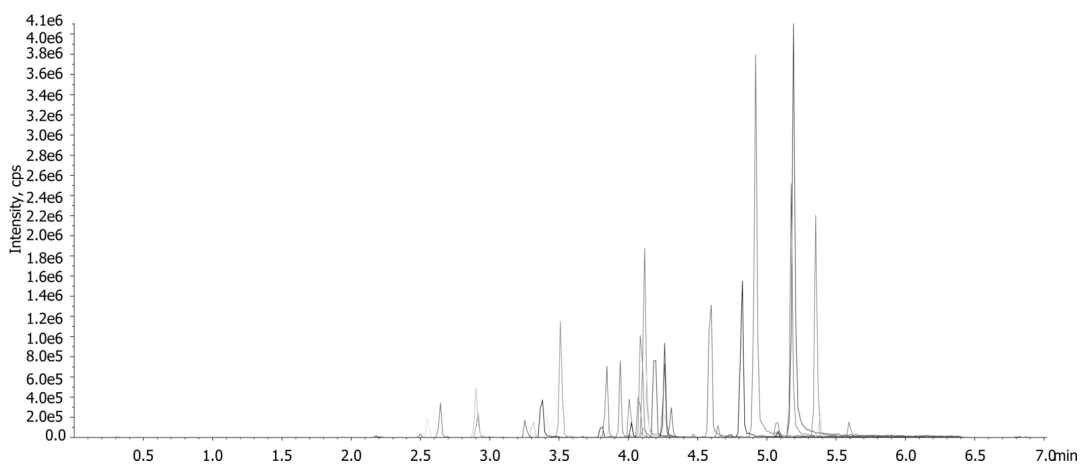

BEH Comps - Pesticide Panel using Luna Omega 1.6um Polar C18 50 x 2.1mm

ANALYTES:

1	Acephate 1	41	Oxamyl 1
2	Acequinocyl 1	42	Paclbutrazol 1
3	Acetamiprid 1	43	Parathion-methyl 1
4	Aldicarb 1	44	Permethrin, cis-1
5	Avermectin B1a 1	45	Permethrin, trans-1
6	Avermectin B1b 1	46	Phosmet 1
7	Azoxystrobin 1	47	Piperonyl butoxide 1
8	Bifenazate 1	48	Prallethrin 1
9	Bifenthrin 1	49	Propiconazole 1
10	Boscalid 1	50	Propoxure 1
11	Carbaryl 1	51	Pyrethrins Cinerin 1
12	Carbofuran 1	52	Pyrethrins Jasmolin 1
13	Chlorantraniliprole 1	53	Pyrethrins Pyrethrin 1
14	Chlorpyrifos 1	54	Pyridaben 1
15	Clofentezine 1	55	Spinosyn A 1
16	Cyfluthrin 1	56	Spinosyn D 1
17	Cypermethrin 1	57	Spiromesifen 1
18	Daminozide 1	58	Spirotetramat 1
19	Diazinon 1	59	Spiroxamine 1
20	Dichlorvos 1	60	Tebuconazole 1
21	Dimethoate 1	61	Thiacloprid 1
22	Ethoprophos 1	62	Thiamethoxam 1
23	Etofenprox 1	63	Trifloxystrobin 1
24	Etoxazole 1	64	Uniconazole 1
25	Fenoxycarb 1		
26	Fenpyroximate 1		
27	Fipronil 1		
28	Fipronil NH4 1		
29	Flonicamid 1		
30	Fludioxinil 1		
31	Hexythiazox 1		
32	Imazalil 1		
33	Imidacloprid 1		
34	Kresoxim-methyl 1		
35	Malathion A 1		
36	Metalaxyl 1		
37	Methiocarb 1		
38	Methomyl 1		
39	Myclobutanil 1		
40	Naled 1		

