

## Chiral Separation of Bupivacaine on Lux 5µm i-Amylose-1 in NP

**Column:** Lux® 5 µm i-Amylose-1, LC Column 250 x 4.6 mm, Ea

**Dimensions:** 250 x 4.6 mm ID

**Order No:** 00G-4762-E0

**Elution Type:** Isocratic

**Eluent A:** 0.1% DEA in Hexane

**Eluent B:** 0.1% DEA in THF

Gradient Profile:	Step No.	Time (min)	Pct A	Pct B
	1	0	85	15


**Flow Rate:** 1 mL/min

**Col. Temp.:** ambient

**Detection:** UV-Vis Abs.-Variable Wave.(UV) @ 254 nm (ambient)  
24503



Products used in this application:

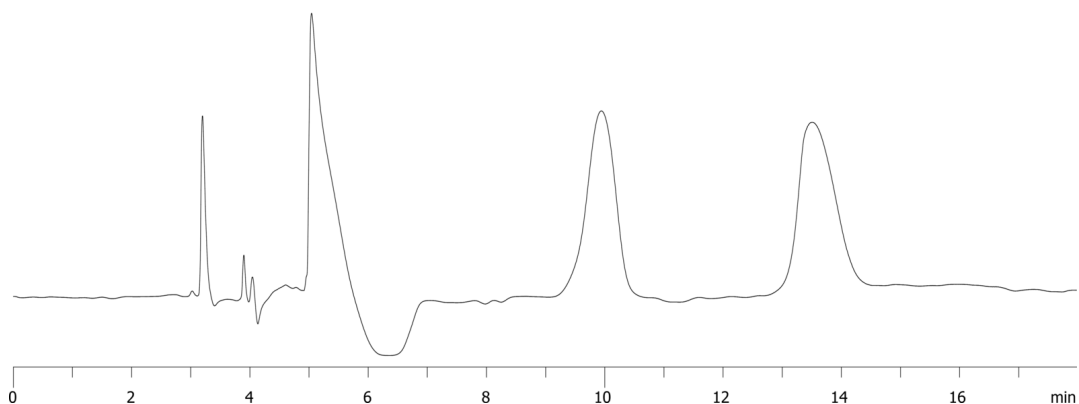

### ANALYTES:

**1** Bupivacaine

Retention Time: 9.941 min

**2** Bupivacaine

Retention Time: 13.497 min

