

## Lux i-Amylose-3 Separation for Econazole

**Column:** Lux® 5 µm i-Amylose-3, LC Column 250 x 4.6 mm, Ea

**Dimensions:** 250 x 4.6 mm ID

**Order No:** 00G-4779-E0

**Elution Type:** Isocratic

**Eluent A:** Water

**Eluent B:** Acetonitrile

**Eluent C:** Diethylamine

Gradient Profile:	Step No.	Time (min)	Pct A	Pct B	Pct C
	1	18	35	65	0.1

**Flow Rate:** 1 mL/min


**Col. Temp.:** 25 °C

**Detection:** LC/UV (DAD, PDA) @ 254.000000000 nm (nanometers) (ambient)

**Analyst Note:** 2 mg/mL in Ethanol



Products used in this application:

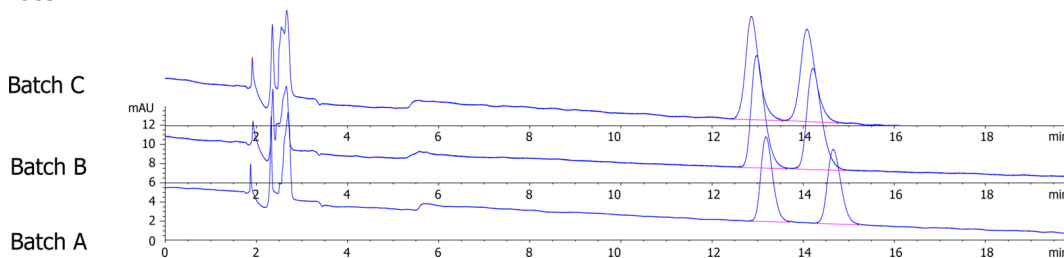

### ANALYTES:

1 Econazole

Retention Time: 13.181 min

2 Econazole

Retention Time: 14.663 min

