

HPLC Application

ID No.: 25372

Lux i-Cellulose-5 Separation for Thalidomide

Column: Lux[®] 5 µm i-Cellulose-5, LC Column 250 x 4.6 mm, Ea

Dimensions: 250 x 4.6 mm ID

Order No: 00G-4756-E0

Elution Type: Isocratic

Eluent A: Acetonitrile

Eluent B: Diethylamine

Gradient Profile:	Step No.	Time (min)	Pct A	Pct B
	1	6	100	0.1


Flow Rate: 1 mL/min

Col. Temp.: 25 °C

Detection: LC/UV (DAD, PDA) @ 254.000000000 nm (nanometers) (ambient)

Analyst Note: 2 mg/mL in Ethanol

Luxi-Cellulose5

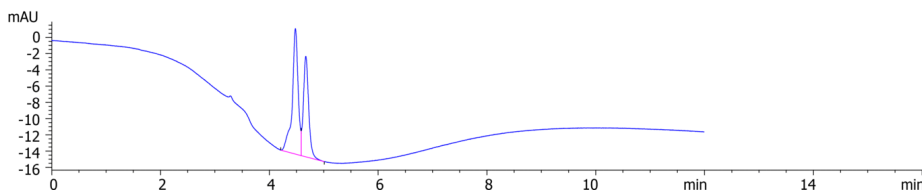

Products used in this application:



ANALYTES:

- 1 Thalidomide
Retention Time: 4.479 min
- 2 Thalidomide
Retention Time: 4.671 min