

# HPLC Application

ID No.: 25653

## Methanol/Acetonitrile (100:0) Impact on Phenyl Selectivity for 2.6 µm Kinetex F5

**Column:** Kinetex<sup>®</sup> F5 2.6µm 100 Å, LC Column 50 x 4.6 mm, Ea

**Dimensions:** 50 x 4.6 mm ID

**Order No:** 00B-4723-E0

**Elution Type:** Gradient

**Eluent A:** 0.1% Formic acid in Water

**Eluent B:** 0.1% Formic acid in Methanol

Gradient Profile:	Step No.	Time (min)	Pct A	Pct B
	1	0	85	15
	2	4	5	95


**Flow Rate:** 800 µL/min

**Col. Temp.:** ambient

**Detection:** LC/MS @ m/z (mass per charge) (ambient)



Products used in this application:

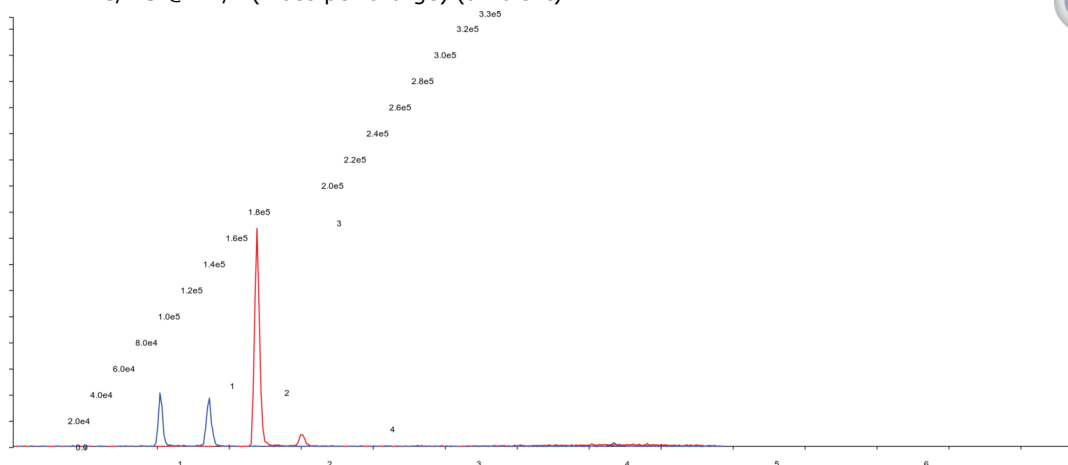

### ANALYTES:

- 1** Morphine  
Retention Time: 1.02 min
- 2** Hydromorphone  
Retention Time: 1.36 min
- 3** Codeine  
Retention Time: 1.69 min
- 4** Hydrocodone  
Retention Time: 2 min