

## Methanol/Acetonitrile (60:40) Impact on Selectivity for Kinetex 2.6 μm C18

**Column:** Kinetex® 2.6 μm C18 100 Å, LC Column 50 x 4.6 mm, Ea

**Dimensions:** 50 x 4.6 mm ID

**Order No:** 00B-4462-E0

**Elution Type:** Gradient

**Eluent A:** 0.1% Formic acid in Water

**Eluent B:** 0.1% Formic acid in Methanol/Acetonitrile (60:40)

Gradient Profile:	Step No.	Time (min)	Pct A	Pct B
	1	0	85	15
	2	4	5	95


**Flow Rate:** 800 μL/min

**Col. Temp.:** ambient

**Detection:** LC/MS @ m/z (mass per charge) (ambient)



Products used in this application:

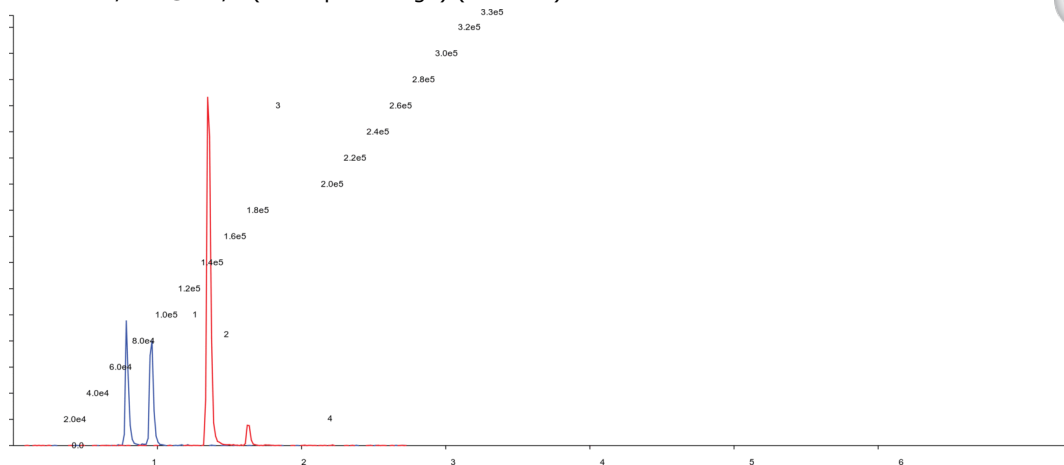

### ANALYTES:

- 1 Morphine**  
Retention Time: 0.79 min
- 2 Hydromorphone**  
Retention Time: 0.96 min
- 3 Codeine**  
Retention Time: 1.35 min
- 4 Hydrocodone**  
Retention Time: 1.63 min

