



GC-MS/MS Separation of Hexachlorodibenzofurans (HxCDF) on ZB-Dioxin GC column with 5m Zebtron Guardian

Column: Zebtron™ ZB-Dioxin w/GUARDIAN™, 5 m. GC Cap. Column 60 m x 0.25 mm x 0.20 µm, E
Phase:
Dimensions: 60 meters x 0.25 mm x 0.2 µm
Order No: 7KG-G045-10-GGA
Oven Profile: to 220°C @ 30°C/min 220°C for 15.5 min to 240°C @ 2°C/min 240°C for 4.5 min to 270°C @ 5°C/min 270°C for 5 min to 330°C @ 15°C/min
Carrier Gas: Constant Flow Helium, 1.056 mL/min
Injection: Splitless :1 1 µL @ 330°C
Detection: GC/MS/MS (340°C)

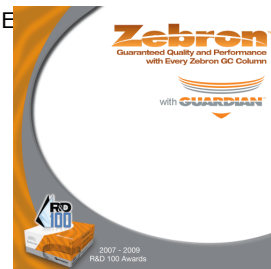

Products used in this application:

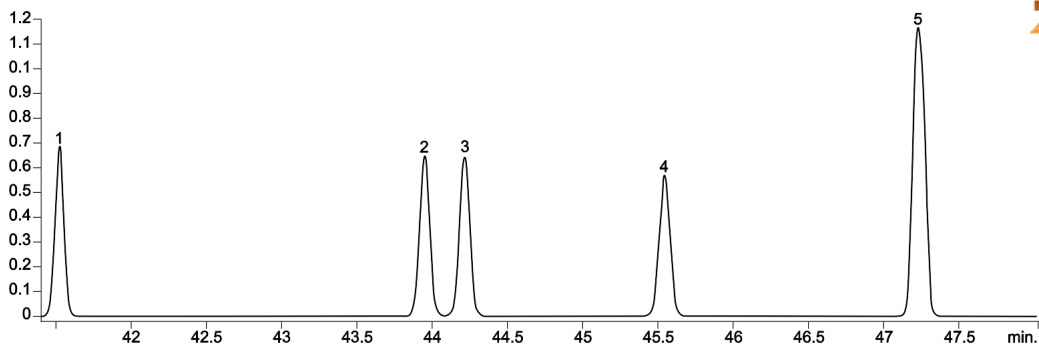

ANALYTES:

- 1** 1,2,3,4,7,8-HxCDF
Retention Time: 41.1 min
- 2** 1,2,3,4,6,8-HxCDF
Retention Time: 43.9 min
- 3** 1,2,3,6,7,8-HxCDF
Retention Time: 44.3 min
- 4** 2,3,4,6,7,8-HxCDF
Retention Time: 45.6 min
- 5** 1,2,3,7,8,9-HxCDF
Retention Time: 47.2 min

