

## GC-MS/MS Separation of Heptachlorodibenzofurans (HpCDF) on ZB-Dioxin GC column with 5m Zebron Guardian

**Column:** Zebron™ ZB-Dioxin w/GUARDIAN™, 5 m. GC Cap. Column 60 m x 0.25 mm x 0.20 µm, E

**Phase:**

**Dimensions:** 60 meters x 0.25 mm x 0.2 µm


**Order No:** 7KG-G045-10-GGA

**Oven Profile:** to 220°C @ 30°C/min 220°C for 15.5 min to 240°C @ 2°C/min 240°C for 4.5 min to 270°C @ 5°C/min 270°C for 5 min to 330°C @ 15°C/min


**Carrier Gas:** Constant Flow Helium, 1.056 mL/min

**Injection:** Splitless :1 1 µL @ 330°C

**Detection:** GC/MS/MS (340°C)

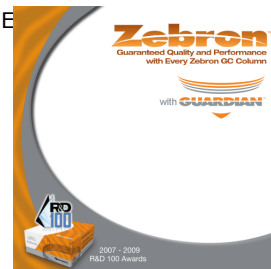

Products used in this application:

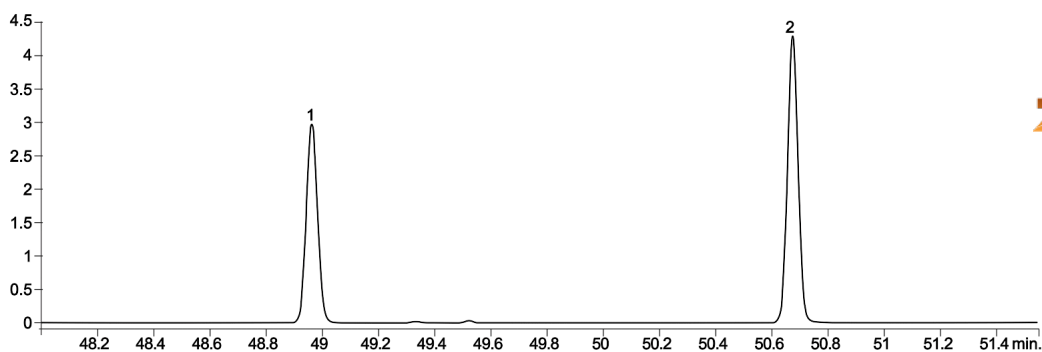

### ANALYTES:

**1** 1,2,3,4,6,7,8-HpCDF

Retention Time: 48.9 min

**2** 1,2,3,4,7,8,9-HpCDF

Retention Time: 50.7 min

