

# HPLC Application

ID No.: 3288

## Pindolol - USP Method

**Column:** Luna® 3 µm C18(2) 100 Å, LC Column 150 x 4.6 mm, Ea

**Dimensions:** 150 x 4.6 mm ID

**Order No:** 00F-4251-E0

**Elution Type:** Isocratic

**Eluent A:** 50mM Sodium acetate, pH 5.0

**Eluent B:** Acetonitrile

<b>Gradient Profile:</b>	<b>Step No.</b>	<b>Time (min)</b>	<b>Pct A</b>	<b>Pct B</b>
<b>1</b>	<b>1</b>	<b>0</b>	<b>65</b>	<b>35</b>

**Flow Rate:** 1 mL/min

**Col. Temp.:** ambient

**Detection:** UV-Vis Abs.-Variable Wave.(UV) @ 219 nm (ambient)


**Analyst Note:** SecurityGuard™ Guard Cartridge System extends column lifetime.

- SecurityGuard Cartridges, C18 4 x 3.0mm, 10/Pk Part No.: AJ0-4287

- Holder Part No.: KJ0-4282



Products used in this application:

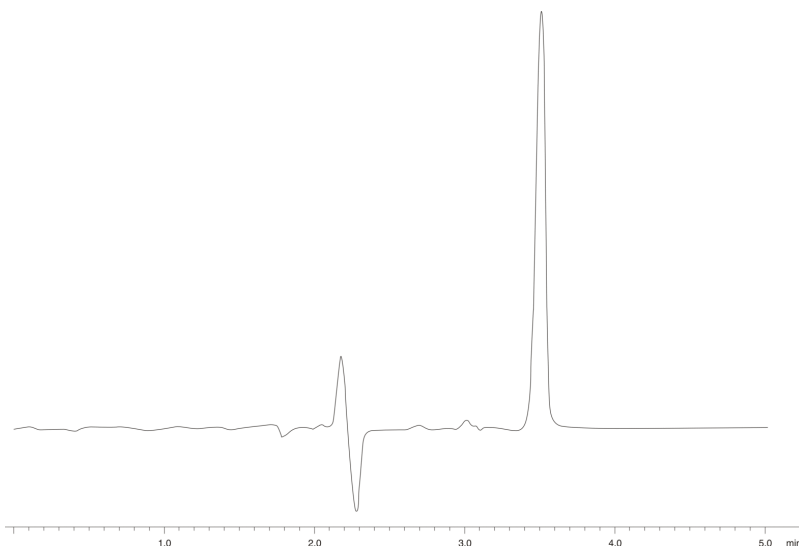

### ANALYTES:

**1** Pindolol

