

# HPLC Application

ID No.: 3309

## Alprazolam - Alternative Method

**Column:** Luna<sup>®</sup> 3  $\mu$ m CN 100 Å, LC Column 100 x 4.6 mm, Ea

**Dimensions:** 100 x 4.6 mm ID

**Order No:** 00D-4254-E0

**Elution Type:** Isocratic

**Eluent A:** Hexane

**Eluent B:** Methylene chloride : Methanol (80:20)

<b>Gradient Profile:</b>	<b>Step No.</b>	<b>Time (min)</b>	<b>Pct A</b>	<b>Pct B</b>
	<b>1</b>	0	75	25

**Flow Rate:** 1 mL/min

**Col. Temp.:** ambient

**Detection:** UV-Vis Abs.-Variable Wave.(UV) @ 254 nm (ambient)


**Analyst Note:** SecurityGuard<sup>™</sup> Guard Cartridge System extends column lifetime.

- SecurityGuard Cartridges, CN 4 x 3.0mm, 10/Pk Part No.: AJ0-4305

- Holder Part No.: KJ0-4282



Products used in this application:

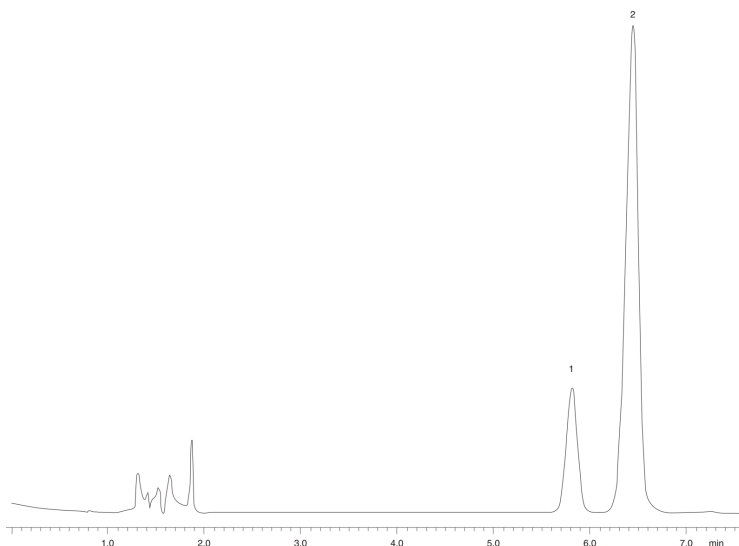

### ANALYTES:

- 1 Triazolam
- 2 Alprazolam

