

## PAH with Surrogate Standards on Envirosep-PP (UV), EPA Method 610


**Column:** N/A**Dimensions:** N/A**Order No:** N/A**Elution Type:** N/A**Eluent A:** Water**Eluent B:** Acetonitrile

Gradient Profile:	Step No.	Time (min)	Pct A	Pct B
	1	0	60	40
	2	25	0	100

**Flow Rate:** 2 mL/min**Col. Temp.:** ambient**Detection:** UV-Vis Abs.-Variable Wave.(UV) @ 254 nm (ambient)**Analyst Note:** SecurityGuard™ Guard Cartridge System extends column lifetime.

- SecurityGuard Cartridges, C18 4 x 3.0mm, 10/Pk Part No.: AJ0-4287

- Holder Part No.: KJ0-4282

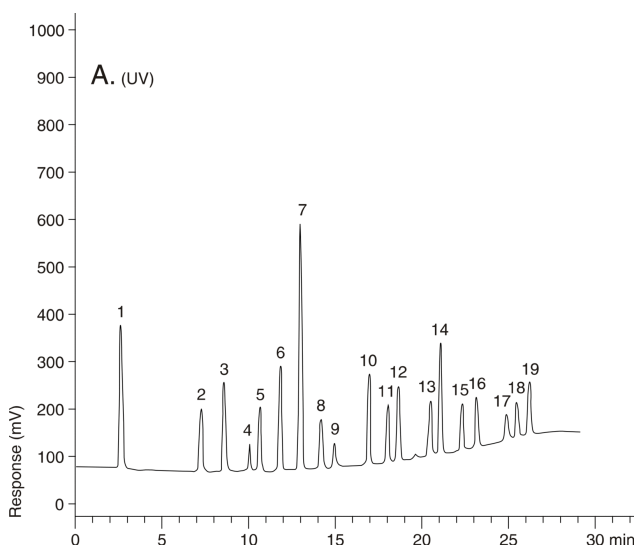

## PAH with Surrogate Standards on Envirosep-PP (UV), EPA Method 610

### ANALYTES:

- 1 Nitrobenzene
- 2 Naphthalene
- 3 Acenaphthylene
- 4 Acenaphthene
- 5 Fluorene
- 6 Phenanthrene
- 7 Anthracene
- 8 Fluoranthene
- 9 Pyrene
- 10 p-Terphenyl
- 11 Benzo[a]anthracene
- 12 Chrysene
- 13 Benzo[e]pyrene
- 14 Benzo[b]fluoranthene
- 15 Benzo[k]fluoranthene
- 16 Benzo[a]pyrene
- 17 Dibenzo[a,h]anthracene
- 18 Benzo[g,h,i]perylene
- 19 Indeno[1,2,3-cd]pyrene

