

## PAH with Surrogate Standards on Envirosep-PP (Flu), EPA Method 610


**Column:** N/A**Dimensions:** N/A**Order No:** N/A**Elution Type:** N/A**Eluent A:** Water**Eluent B:** Acetonitrile

Gradient Profile:	Step No.	Time (min)	Pct A	Pct B
	1	0	60	40
	2	25	0	100

**Flow Rate:** 2 mL/min**Col. Temp.:** ambient**Detection:** Fluorescence (FLUOR) @ 240 nm (ambient)**Analyst Note:** SecurityGuard™ Guard Cartridge System extends column lifetime.

- SecurityGuard Cartridges, C18 4 x 3.0mm, 10/Pk Part No.: AJ0-4287

- Holder Part No.: KJ0-4282

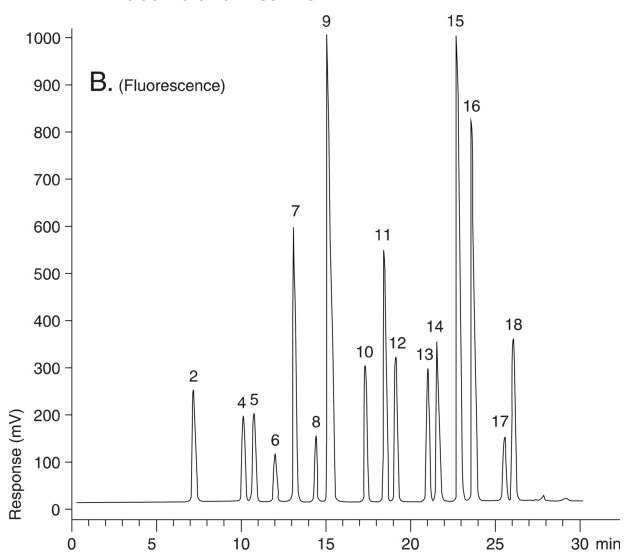

## PAH with Surrogate Standards on Envirosep-PP (Flu), EPA Method 610

### ANALYTES:

- 1** Nitrobenzene
- 2** Naphthalene
- 3** Acenaphthylene
- 4** Acenaphthene
- 5** Fluorene
- 6** Phenanthrene
- 7** Anthracene
- 8** Fluoranthene
- 9** Pyrene
- 10** p-Terphenyl
- 11** Benzo[a]anthracene
- 12** Chrysene
- 13** Benzo[e]pyrene
- 14** Benzo[b]fluoranthene
- 15** Benzo[k]fluoranthene
- 16** Benzo[a]pyrene
- 17** Dibenzo[a,h]anthracene
- 18** Benzo[g,h,i]perylene
- 19** Indeno[1,2,3-cd]pyrene

