

## Peptide Hormones

**Column:** Jupiter® 5 µm C4 300 Å, LC Column 250 x 4.6 mm, Ea

**Dimensions:** 250 x 4.6 mm ID

**Order No:** 00G-4167-E0

**Elution Type:** Gradient

**Eluent A:** 0.1% TFA in water

**Eluent B:** 0.1% TFA in acetonitrile/water (90:10)

Gradient Profile:	Step No.	Time (min)	Pct A	Pct B
	1	0	90	10
	2	8	74	26

**Flow Rate:** 1 mL/min

**Col. Temp.:** ambient

**Detection:** UV-Vis Abs.-Variable Wave.(UV) @ nm (ambient)


**Analyst Note:** SecurityGuard™ Guard Cartridge System extends column lifetime.

- SecurityGuard Cartridges, Widedpore C4 4 x 3.0mm, 10/Pk Part No.: AJ0-4330

- Holder Part No.: KJ0-4282



Products used in this application:

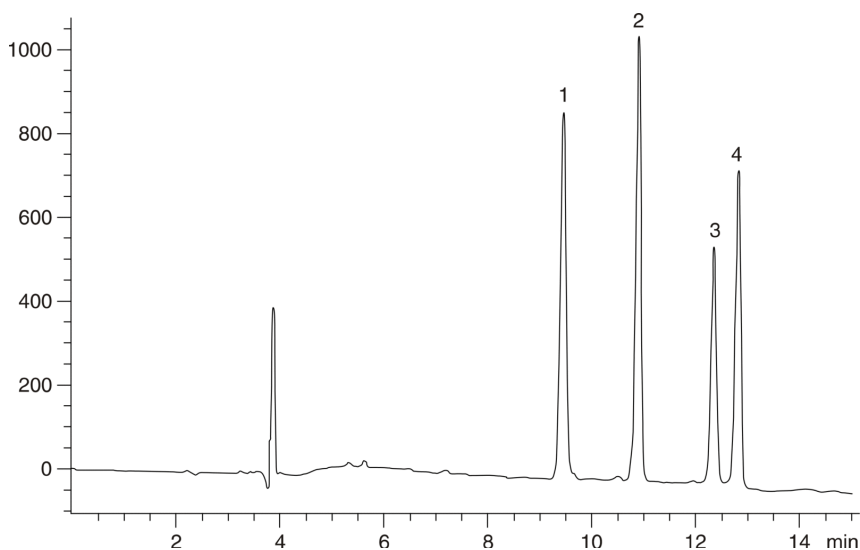

### ANALYTES:

- 1 [Arg8]-Vasotocin
- 2 [Arg8]-Vasopressin
- 3 Isotocin
- 4 Oxytocin

